# NESLAB Merlin Recirculating Chillers

Thermo Manual P/N U00753  
Rev. 02/24/04

**Installation  
Operation  
Basic Maintenance**

No part of this publication may be reproduced or reprinted except for personal use without the prior permission of Thermo Electron.

**Thermo**  
ELECTRON CORPORATION

# Table of Contents

## Preface

Compliance .....	3
After-sale support.....	3
Feedback .....	4
Warranty .....	4
NES-care Extended Warranty Contract .....	4
Unpacking .....	4

## Section I Safety

Warnings .....	5
Additional Warnings .....	6

## Section II General Information

Quick Reference Operating Procedures .....	7
Description .....	16
Specifications .....	16

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Site .....	21
Electrical Requirements .....	22
Plumbing Requirements .....	24
Fluids .....	25
Filling Requirements .....	25
Reservoir Tank Isolation Valves/Open Tank Circulation .....	26
Controller .....	27
Start Up/Shut Down .....	29
Setup/Tuning Loop .....	31

## Section IV Options

Communication Connector .....	33
Remote Start/Stop .....	34
External Temperature Sensor .....	35
Reservoir Level Warning/Fault .....	36
Low Pump Flow Warning/Fault.....	36
Additional Accessories .....	36

## Section V Maintenance

Error Codes .....	37
Displaying Software Version Number .....	38
Controller PID Values .....	38
Cleaning .....	39
Algae .....	39
Leaks .....	39
Phase Rotation Interlock .....	40
PD Pressure Relief Valve .....	40
Pump Strainer .....	41
PD Pump Motor Lubrication .....	42
Storage .....	42
BOM Decoder .....	42
Temperature Sensor Calibration .....	43

## Section VI Troubleshooting

Checklist .....	45
Service Assistance .....	46

## Appendix A Water Quality Standards and Recommendations

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## Appendix C RS-485 NC Serial Communications Protocol

Blank page.

## Preface

### Compliance

#### 60 Hertz

Products tested and found to be in compliance with the requirements defined in CSA-C22.2 No 1010-1 and UL-3101-1 can be identified by the C(UL)US label on the unit's serial tag.

#### 50 Hertz

Products tested and found to be in compliance with the requirements defined in the EMC standards defined by 89/336/EEC as well as Low Voltage Directive (LVD) 73/23/EEC can be identified by the CE Mark on the rear of the unit. The testing has demonstrated compliance with the following directives:

LVD, 73/23/EEC	EN 61010-1:1993
EMC, 89/336/EEC	EN 55011, Class A Verification
	EN 50082-1:1992
	IEC 61000-4-2:1999
	IEC 61000-4-3:1998
	IEC 61000-4-4:1995
	IEC 61000-4-5:1995
	IEC 61000-4-6:1996

For any additional information, refer to the Declaration of Conformity that shipped with the unit.

### After-sale Support

Thermo Electron Corporation is committed to service both during and after the sale. If you have questions concerning operation, contact our Sales Department. If your unit fails to operate properly, or if you have questions concerning spare parts or Service Contracts, contact our Service Department.

Before calling, *please* obtain the following information:

- unit BOM number \_\_\_\_\_
- unit serial number \_\_\_\_\_
- unit software version (page 38) \_\_\_\_\_
- voltage of power source \_\_\_\_\_

The unit's BOM and serial number are located on the name plate label on the rear of the unit. See page 42 for instructions on how to decode your unit's BOM number.

Thermo Electron			
Newington, NH 03801 U.S.A.			
(800)258-0830 / (603)436-9444			
<b>BOM#:</b>	XXXXXXXXXXXX		
<b>S/N:</b>	XXXXXXXXXX		
XXX VOLT	XX HZ	X PH	X.X AMP
R134A	X.XOZ	HIGH XXX PSI	LOW XXXPSIG

## **Feedback**

We appreciate any feedback you can give us on this manual. Please e-mail us at [neslabmanuals@thermo.com](mailto:neslabmanuals@thermo.com). Be sure to include the manual part number and the revision date listed on the front cover.

## **Warranty**

Units have a warranty against defective parts and workmanship for 24 months from date of shipment. See back page for more details.

## ***NES-care* Extended Warranty Contract**

- Extend parts and labor coverage for an additional year.
- Worry-free operation.
- Control service costs.
- Eliminate the need to generate repair orders.
- No unexpected repair costs.

Other contract options are available. Please contact Thermo for more information.

## **Unpacking**

Retain all cartons and packing material until the unit is operated and found to be in good condition. If the unit shows external or internal damage contact the transportation company and file a damage claim. Under ICC regulations, this is your responsibility.

## Section I Safety

### Warnings

Warnings are posted throughout the manual. These warnings are designated by an exclamation mark inside an equilateral triangle and text highlighted in bold. Read and follow these important instructions. Failure to observe these instructions can result in permanent damage to the unit, significant property damage, or personal injury or death.

The lightning flash with arrow symbol, within an equilateral triangle, is intended to alert the user to the presence of non-insulated "dangerous voltage" within the unit's enclosure. The voltage may be of significant magnitude to constitute a risk of electrical shock.

Make sure you read and understand all instructions and safety precautions listed in this manual before installing or operating your unit. If you have any questions concerning the operation of your unit or the information in this manual, please contact our Sales Department (see After-sale Support).

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**The unit construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection may not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided.**

**Never connect the inlet or outlet fitting to your building water supply or any water pressure source.**

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**Never use flammable or corrosive fluids with this unit. Highly distilled and deionized water may be aggressive and cause material corrosion. Please contact Thermo before subjecting this unit to prolonged exposure to highly distilled or deionized water.**

**Do not use automotive antifreeze. Commercial antifreeze contains silicates that can damage the pump seals. Use of automotive antifreeze will void the manufacturer's warranty.**

## **Additional Warnings**

In addition to the specific warnings listed on the previous page the following general warnings apply to your unit:

**Performance of installation, operation, or maintenance procedures other than those described in this manual may result in a hazardous situation and may void the manufacturer's warranty.**

**Transport the unit with care. Sudden jolts or drops can damage the refrigeration lines.**

**If the unit is to be transported and/or stored it needs to be drained. If the unit is to be transported and/or stored in cold temperatures it needs to be drained and then flushed with a 50/50 glycol/water mixture.**

**Observe all warning labels.**

**Never remove warning labels.**

**Never operate damaged or leaking equipment.**

**Never operate the unit without cooling fluid in the reservoir.**

**Always turn off the unit and disconnect the power cord from the power source before performing any service or maintenance procedures, or before moving the unit.**

**Never operate equipment with damaged power cords.**

**Refer service and repairs to a qualified technician.**

## Section II General Information


Indicates **refrigeration system status**. See page 27.

**Optional.** Flashes to indicate the **external sensor** is enabled. See page 35.

Indicates the controller is displaying the **setpoint**. Press the arrow keys to change the value. See page 30.

Indicates the controller is displaying the unit's **low temperature alarm setting**. Press the arrow keys to change the value. See page 30.

Press to **sequence** through the four indicators and accept the displayed value. See page 27.  
**NOTE:** None of the four indicators are visible unless this key is pressed.

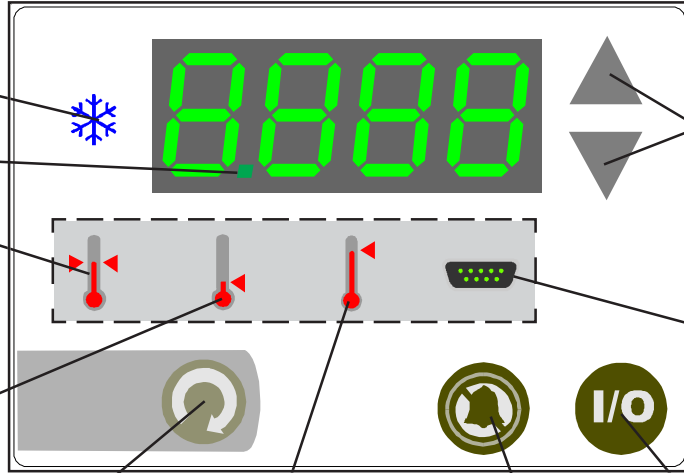Press to **change** the displayed value. See page 27.


**Optional.** Indicates the unit is operating in the **serial communication** or remote start/stop mode. See page 31.

Press to **start/stop** the unit. See page 27.

Indicates the controller is displaying the **high temperature alarm setting**. Press the arrow keys to change the value. See page 30.

Press to **mute** the alarm. See page 27.

Reservoir Access Panel. See page 24.

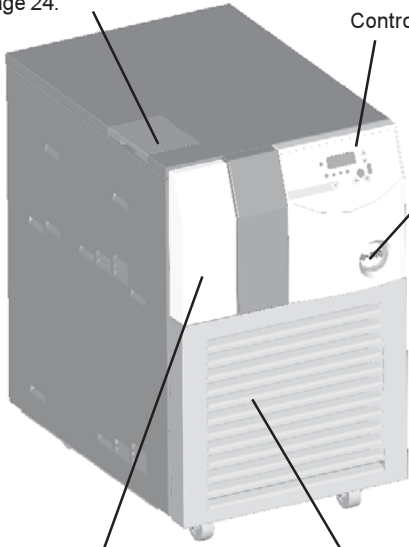Controller. See page 27.

Pressure Gauge. See page 29.

Indicates a **Merlin Operating Tip**.

Reservoir. See page 25.

Condenser Access Panel. See page 39.

Circuit Breaker. See page 27.

**Optional**  
Communications Connector. See page 33.

Plumbing Return Connection. See page 24.

Plumbing Supply Connection. See page 24.

PSI/KPa PD Pump Pressure Relief Valve. See page 40.

**Optional**  
Remote Start/Stop Connector. See page 34.  
External Temperature Sensor. See page 35

Drain. See page 24.

## Quick Reference Operating Procedures

### Installation

The unit has an air-cooled refrigeration system. Air is drawn in the front of the unit and discharged through rear and sides. Position the unit so the intake and discharge are not impeded. Inadequate ventilation will reduce cooling capacity and, in extreme cases, result in compressor failure.

Excessively dusty areas should be avoided and a periodic cleaning schedule should be instituted. For proper operation, the unit needs to pull substantial amounts of air through a condenser. A build up of dust or debris on the fins of the condenser will lead to a loss of cooling capacity.

The unit will retain its full rated capacity in ambient temperatures up to approximately +77°F (+25°C).

Make sure the voltage of the power source meets the specified voltage,  $\pm 10\%$ .

The plumbing connections are located on the rear of the unit and are labelled  and . These connections are 1/2 inch FPT (3/4 inch FPT for units with CP-55 pumps). Remove the plastic protective plugs from both plumbing connections. Connect the  fitting to the inlet of your application. Connect the  fitting to the outlet of your application.

To fill the reservoir, remove the reservoir access panel and locate the reservoir cap. Remove the cap and fill the reservoir with clean cooling fluid.

For fluid recommendations refer to page 25.

### Operation

Before starting the unit, double-check all electrical and plumbing connections. Make sure the circulating system is properly filled with cooling fluid. Place the circuit breaker located on the rear of the unit to the up position.

To start the unit, press . To turn the unit off press  again.

The  LED indicates the status of the refrigeration system. It illuminates to indicate the refrigeration system is removing heat from the cooling fluid. As the operating temperature approaches the setpoint, the LED will flash.

### Temperature Adjustment

To display the temperature setpoint, press  on the controller. The  indicator will illuminate and the display will flash the current setpoint value. To adjust the temperature setpoint, press the arrow keys until the desired temperature setpoint is indicated. Press  again to confirm the change. The display will rapidly flash the new value for a short time and then return to the recirculating fluid temperature.

### Periodic Maintenance

Periodically inspect the reservoir fluid. If cleaning is necessary, flush the reservoir with a cleaning fluid compatible with the circulating system and the cooling fluid.

The cooling fluid should be replaced periodically. Frequency depends on the operating environment and amount of usage.

Before changing the cooling fluid ensure it is at a safe temperature.

Periodic vacuuming of the condenser fins is necessary. The frequency of cleaning depends on the operating environment. We recommend a monthly visual inspection of the condenser after initial installation. After several months, the cleaning frequency will be established.

Units with PD pumps have a strainer. If debris is in the system, the strainer will prevent the material from being drawn into the pump and damaging the pump vanes. A clogged strainer will also cause increased pump discharge pressure.

After initial installation, the strainer may become clogged.  **Clean the strainer after the first week of installation.** After this first cleaning, we recommend a monthly visual inspection. After several months, the cleaning frequency will be established. Before cleaning the strainer, disconnect the power cord from the power source and drain the unit.

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## Description

The NESLAB Merlin Recirculating Chiller is designed to provide a continuous supply of cooling fluid at a constant temperature and volume. The unit consists of an air-cooled refrigeration system, plate heat exchanger, recirculating pump, polyethylene reservoir, and a microprocessor controller.





Throughout the manual, you will be asked to consult the unit's name plate label for specific information. The label is located on the rear of the unit. See page 3 for an example.

## Specifications

	M 25	M 33	M 75	M 100	M 150
Temperature Range <sup>1,2</sup>	+5°C to +35°C				
Temperature Stability <sup>1,3,4</sup>	±0.1°C		±0.15°C		
Reservoir Volume					
Gallons	0.5				
Liters	1.8				
Refrigerant	R134a			R404a	
Cooling Capacity (Watts) <sup>1,4</sup> Standard Temp Units	<div style="display: flex; flex-direction: column; align-items: center;"> <div style="display: flex; justify-content: space-between; width: 100%;"> <span>Fluid - Water</span> <div style="font-size: small;">                     A = M 75 60Hz                      B = M 75 50Hz                      C = M 33 60Hz                      D = M 33 50Hz                      E = M 25 60Hz                      F = M 25, 50Hz                 </div> </div> <div style="display: flex; justify-content: space-between; width: 100%;"> <span>Fluid - Water</span> <div style="font-size: small;">                     A = M 150 60Hz                      B = M 150 50Hz                      C = M 100 50 Hz, 60Hz                 </div> </div> </div>				

1. Specifications subject to change.
2. A temperature range of -15°C to +35°C is available.
3. At the chiller. Display resolution 1.0°C. 0.1°C display resolution is available using the controller's Setup/Tuning Loop, see page 31.
4. Circulating at 20°C ambient, pump unloaded. Cooling capacity will vary depending on fluid temperature, ambient temperature, and cooling fluid.

**Cooling Capacity (Watts)<sup>1,2</sup>**  
**Low Temp Units**

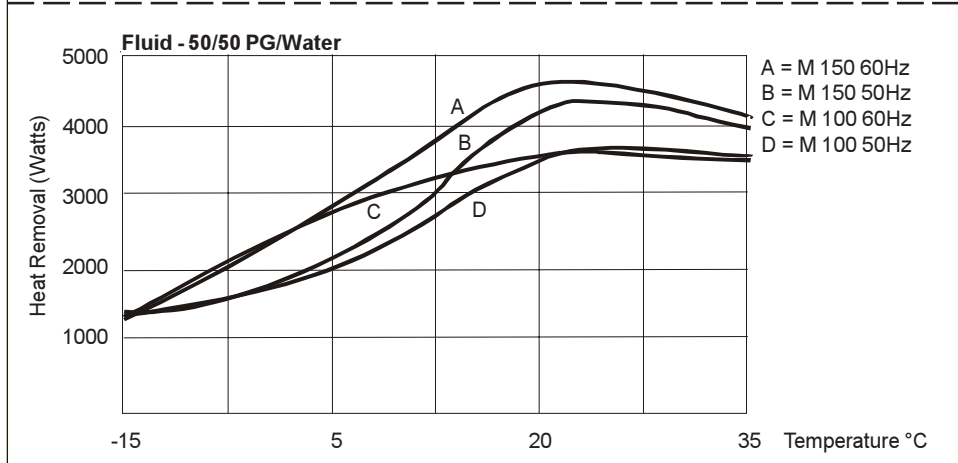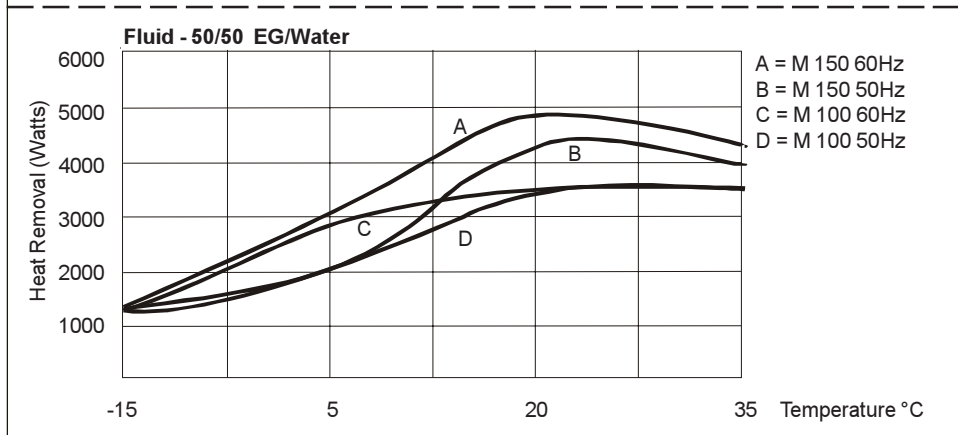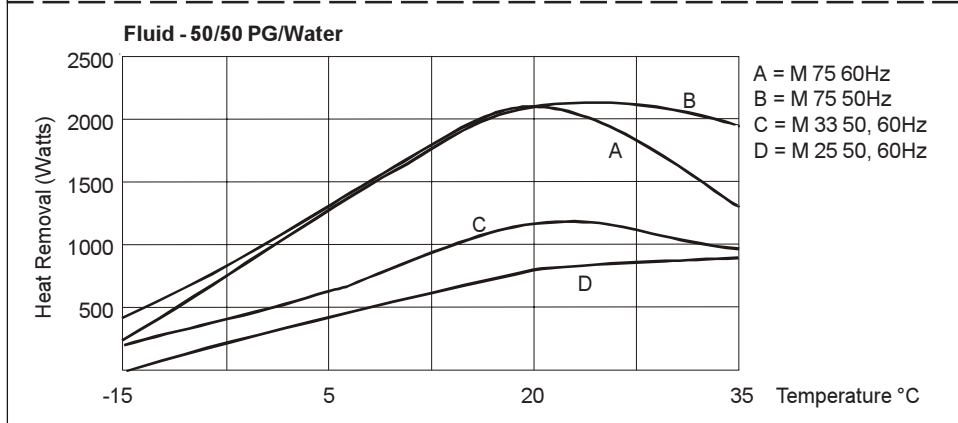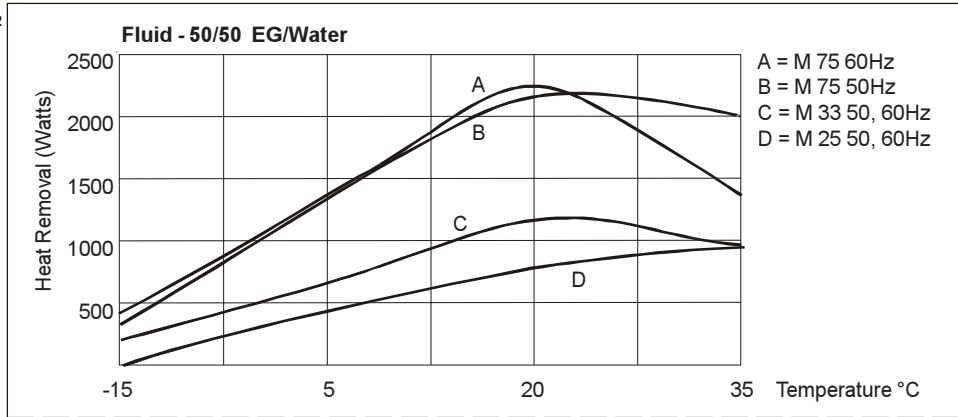1. Specifications subject to change.
2. Circulating at 20°C ambient, pump unloaded. Cooling capacity will vary depending on fluid temperature, ambient temperature, and cooling fluid.

## Specifications

**Intake Requirements (ft<sup>3</sup>/min)**  
(60 Hertz units)

M 25	M 33	M 75	M 100	M 150
178	178	385	755	755
<b>Shipping Weight</b>				
(Pounds)				
141	141	200	287 <sup>1</sup>	287 <sup>1</sup>
(Kilograms)				
64	64	91	130	130

1. Units with PD 2 pumps. Units with CP 55 pumps weigh 310 pounds (141 kilograms).

## Pump Options

### Standard temperature range units

Unit	Standard	Optional
M 25	MD 30 or PD 1	None
M 33	MD 30 or PD 1	PD 2
M 75	PD 1	PD 2
M 100(208/230V 60Hz)	PD 2	CP 55
M 100(230V 50Hz)	PD 2	None
M 100(400V 50Hz)	PD 2	CP 55
M 150(208/230V 60Hz)	PD 2	CP 55
M 150(230V 50Hz)	PD 2	None
M 150(400V 50Hz)	PD 2	CP 55




### Low temperature range units

Unit	Standard	Optional
M 25	PD 1	None
M 33	PD 1	PD 2
M 75	PD 2	None
M 100(208/230V 60Hz)	PD 2	CP 55
M 100(230V 50Hz)	PD 2	None
M 100(400V 50Hz)	PD 2	CP 55
M 150(208/230V 60Hz)	PD 2	CP 55
M 150(230V 50Hz)	PD 2	None
M 150(400V 50Hz)	PD 2	CP 55


**NOTE:** PD pumps should not be used with applications that are sensitive to high pressures.

# Pump Capacity

Description

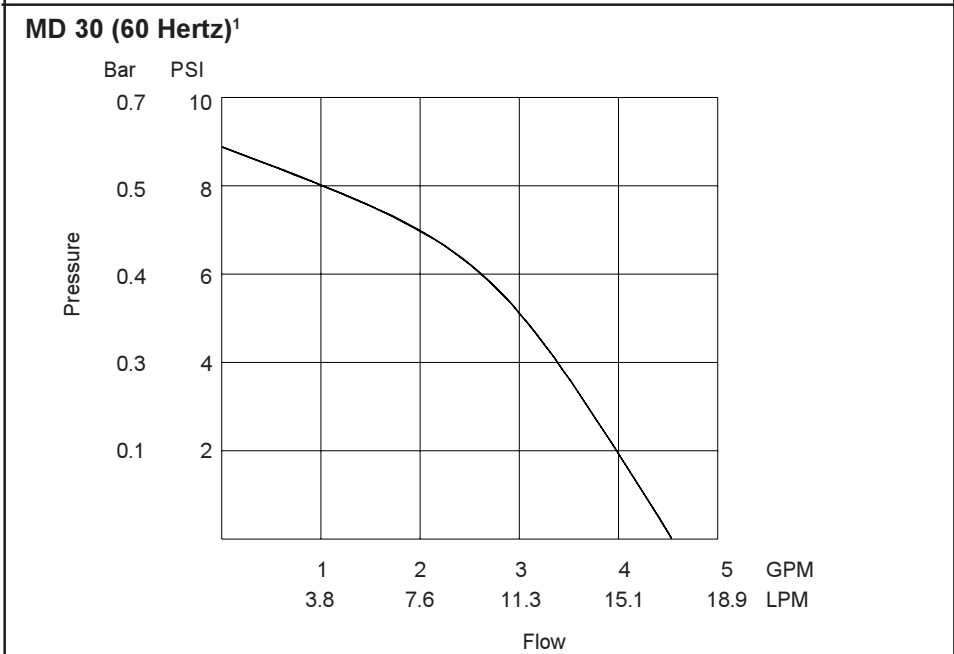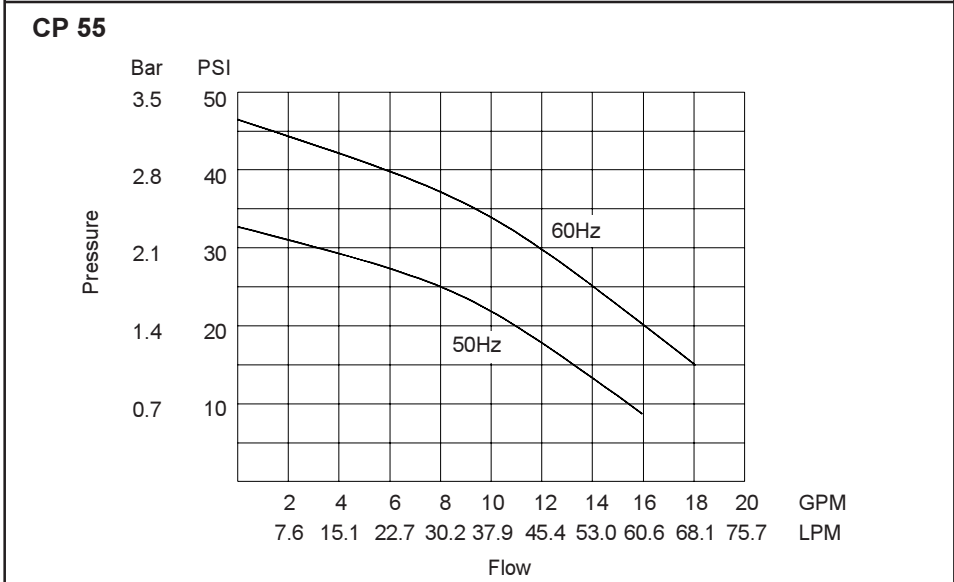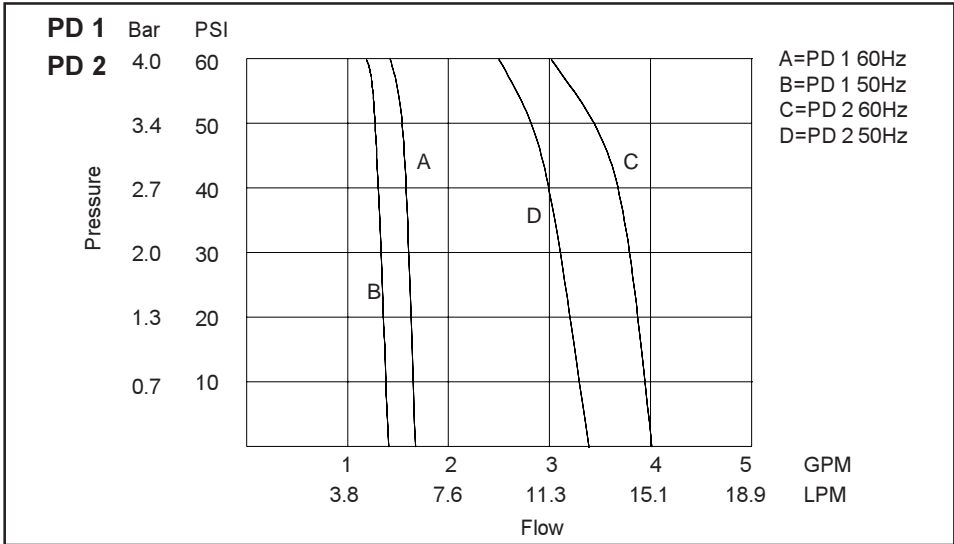1. The MD 30 is not available for 50 Hertz units.

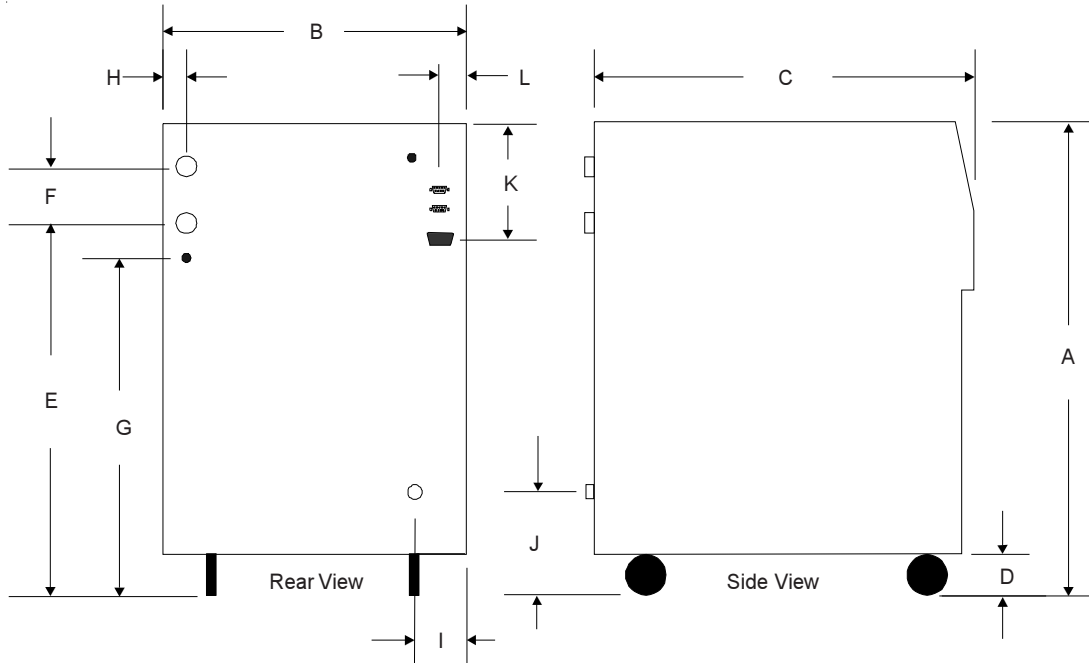Unit Dimensions	M 25	M 33	M 75	M 100	M 150
Dimension A <sup>1</sup>	23 1/2	23 1/2	26 1/8	30 3/8	30 3/8
Dimension B <sup>1</sup>	12 5/8	12 5/8	16 1/4	21 1/4	21 1/4
Dimension C <sup>1</sup>	20 7/8	20 7/8	24	29 1/4	29 1/4
Dimension D <sup>2</sup>	2	2	2	3	3
Dimension E <sup>3</sup>	18 1/8	18 1/8	20 3/4	25	25
Dimension F <sup>4</sup>	3	3	3	3	3
Dimension G <sup>5</sup>	16 1/4	16 1/4	18 7/8	23	23
Dimension H <sup>6</sup>	1 3/8	1 3/8	1 3/8	1 3/8	1 3/8
Dimension I <sup>7</sup>	2 7/8	2 7/8	3 5/8	4 3/4	4 3/4
Dimension J	5	5	7	7 1/2	7 1/2
Dimension K	7 1/8	7 1/8	7 1/8	7 1/8	7 1/8
Dimension L	1 3/4	1 3/4	1 3/4	1 3/4	1 3/4
Crate Dimensions (H x W x D) <sup>8</sup>	36 x 18 1/2 x 26	36 x 18 1/2 x 26	39 x 22 1/4 x 29 1/4	42 x 27 x 35 1/4	42 x 27 x 35 1/4

1. Dimension A is the height of the unit. Units with CP-55 pumps are 6" higher. The unit width and depth (dimensions B and C) are the case dimensions, add approximately 1/2 inch to include the plumbing connections.
2. Dimension D is the distance from the floor to the bottom of the unit case (height of the casters).
3. Dimension E is the distance from the floor to the center of the outlet connection.
4. Dimension F is the distance between the center of the outlet and inlet connections.
5. Dimension G is the distance from the floor to the center of the pressure relief valve (units with PD pumps).
6. Dimension H is the distance from the unit's side to the center of the outlet/inlet connections, and the pressure relief valve. 2" from the unit's side to the center of the outlet/inlet connections for M-100/150 units with CP-55 pumps.
7. Dimension I is the distance from the unit's side to the center of the drain connection. 5 1/4" for M-100/150 units with CP-55 pumps.
8. Dimension J is the distance from the floor to the center of the drain connection. 7 1/2" for M-100/150 units with CP-55 pumps.
9. Dimension K is the distance from the top of the unit to the center of the power connection.
10. Dimension L is the distance from the unit's side to the center of the power connection. NA for M-100/150 units rated 400V/50Hz.

## Section III Installation and Operation

### Site

The unit is designed for indoor use only. It should be located in a clean industrial environment where ambient temperatures are inside the range of 10°C to 35°C (50°F to 94°F) and relative humidity between 10% to 80%, noncondensing.

Operating altitude is sea level to 6000 feet (1830 meters). The unit's over-voltage category is II, its pollution degree is 2.


**Never place the unit in a location where excessive heat, moisture, or corrosive materials are present.**

The unit has an air-cooled refrigeration system. Air is drawn through the front of the unit and discharged through rear and side panels. The unit must be positioned so the intake and discharge are not impeded. A minimum clearance of 3 feet (1 meter) on all vented sides is necessary for adequate ventilation. Inadequate ventilation will cause a reduction in cooling capacity and, in extreme cases, compressor failure.

Excessively dusty areas should be avoided and a periodic cleaning schedule should be instituted (see Section IV, Maintenance and Service). Optional air filters are available, contact our Service Department. See Preface, After-Sale Support.

The unit retains its full rated capacity in ambient temperatures up to approximately 25°C (77°F). Reduce the cooling capacity 1% for every 0.5°C (1°F) above 25°C (77°F), up to a maximum ambient temperature of 35°C (94°F).

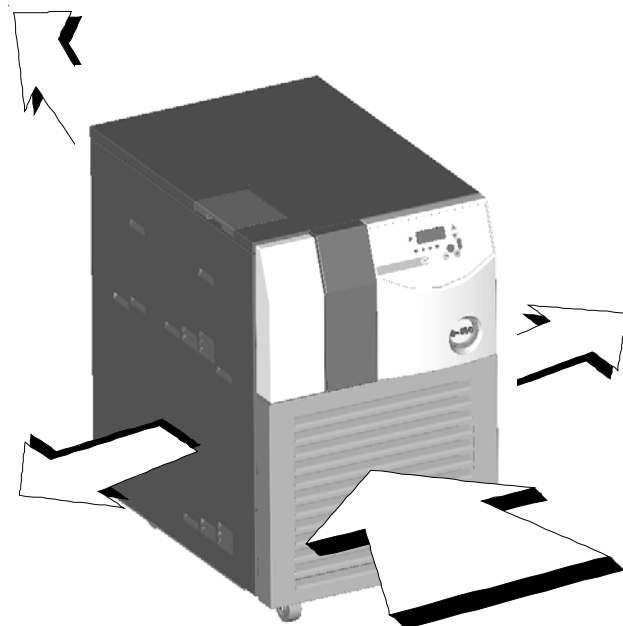Direction of Air Flow

## Electrical Requirements

The unit construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection may not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided.

Available power options:

Standard and Low-Temperature 60 Hertz Units				
Unit	Pump	Voltage	Current Value	Line Cord Plug
M 25	MD 30	115	8.1	5-15P
	PD 1	115	9.8	5-15P
M 33	MD 30	115	11.5	5-20P
	PD 1	115	13.2	5-20P
	PD 2	115	15.5	5-20P
	PD 1	208/230	6.5	6-15P
	PD 2	208/230	7.8	6-15P
M 75	PD 1	208/230	8.8	6-15P
	PD 2	208/230	10.1	6-15P
M 100	PD 2	208/230	13.1	L6-30P
	CP 55	208/230	16.1	L6-30P
M 150	PD 2	208/230	15.4	L6-30P
	CP 55	208/230	18.4	L6-30P

**NOTE:** All 60 Hertz units are UL Listed and Certified to Canadian Standards.

Amperage based on extreme operating conditions. Operating in a 20°C ambient and at nominal voltage will result in a lower amperage draw.

Refer to the name plate label on the rear of the unit for specific electrical requirements.

**NOTE:** Three-phase units are not supplied with a line cord. Instructions for wiring three-phase units are located behind the power box cover. To access the power box remove the unit's top and left-side panel. Remove the screws securing the cover.

Standard and Low-Temperature 50 Hertz Units				
Unit	Pump	Voltage	Current Value	Line Cord Plug
M 25	PD 1	230	5.1	Country Specific
M 33	PD 1	230	6.3	Country Specific
	PD 2	230	7.6	Country Specific
M 75	PD 1	230	8.8	Country Specific
	PD 2	230	10.1	Country Specific
M 100	PD 2	230	13.5	Country Specific
	PD 2	400/3Ø/5W	7.7	Country Specific
	CP 55	400/3Ø/5W	5.5	Country Specific
M 150	PD 2	230	15.6	Country Specific
	PD 2	400/3Ø/5W	7.7	Country Specific
	CP 55	400/3Ø/5W	5.5	Country Specific

Amperage based on extreme operating conditions. Operating in a 20°C ambient and at nominal voltage will result in a lower amperage draw.

Refer to the name plate label on the rear of the unit for specific electrical requirements.

**NOTE:** Three-phase units are not supplied with a line cord. Instructions for wiring three-phase units are located behind the power box cover. To access the power box remove the unit's top and left-side panel. Remove the screws securing the cover.

## Plumbing Requirements


The plumbing connections are located on the rear of the unit and are labeled  (supply) and  (return). The connections are ½ inch Female Pipe Thread. For units with CP-55 pumps, the connections are ¾ inch Female Pipe Thread. Units with ½ inch fittings are supplied with ⅜ inch and ½ inch barbed adapters, units with ¾ inch fittings are supplied with ½ inch and ¾ inch barbed adapters.

Remove the plastic protective plugs from both plumbing connections. Install the barbed adapters to these connections.

Connect  to the hose feeding the inlet of your application. Connect  to the hose from the outlet of your application. Clamp all connections.

**Never connect the fittings to your building water supply or any water pressure source.**

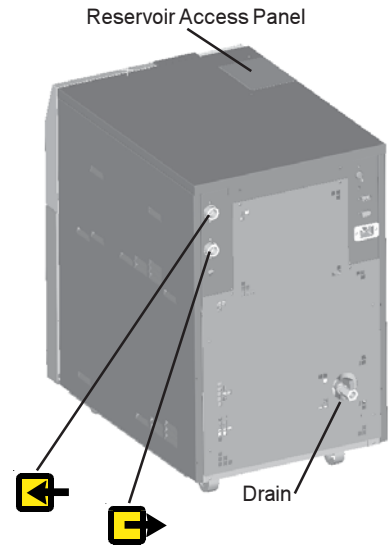**NOTE:** On units equipped with PD pumps, ensure your plumbing is rated to withstand 110 psi at the highest operating temperature.

**On units with MD pumps never completely restrict flow to your application. Dead-heading the pump will damage its coupling and require pump replacement.**

It is important to keep the distance between the unit and the instrument being cooled as short as possible. Tubing should be straight and without bends. If diameter reductions must be made, they should be made at the inlet and outlet of your application, not at the chiller.

All units have a ball-valve reservoir drain located on the rear of the unit. Add the desired length of tubing to the drain valve. To help prevent any spillage from your application, consider installing external isolation valves on the supply and return lines.

## Fluids

Never use flammable or corrosive fluids with this unit. Do not use automotive antifreeze. Commercial antifreeze contains silicates that can damage the pump seals. Use of automotive antifreeze will void the manufacturer's warranty.

Acceptable fluids and their normal operating temperature ranges are:

Filtered/Single Distilled water, +5°C to +35°C (+10°C to +35°C for low-temperature units)

50/50 Uninhibited Ethylene Glycol/Water, -15°C to +35°C

50/50 Inhibited Ethylene Glycol/Water, -15°C to +35°C

50/50 Uninhibited Propylene Glycol/Water, -15°C to +35°C

50/50 Inhibited Propylene Glycol/Water, -15°C to +35°C

Deionized water (1 - 3 megohm, compensated), +5°C to +35°C (+10°C to +35°C for low-temperature units)

Ethylene glycol (EG) is hygroscopic, it will absorb water from its environment. This can affect the freezing point and boiling point of the fluid over time. This may result in system failure.

When using EG/water or PG/water, top-off with EG/water or PG/water. Do not top-off with plain water. Topping-off with plain water can severely effect the freezing point and boiling point of the fluid. This may result in system failure.


Do not use a Deionization (DI) filter with Inhibited EG or Inhibited PG. A DI filter will remove inhibitors from the solution rendering the fluid ineffective against corrosion protection. Also, inhibitors increase fluid conductivity.

Refer to Appendix A for additional information.

## Filling Requirements

The polyethylene reservoir is translucent for easy fluid level monitoring. Remove the reservoir access panel, see previous page. Locate and remove the reservoir cap. Fill the reservoir with clean cooling fluid.


Since the reservoir capacity may be small compared to your application, have extra cooling fluid on hand to keep the system topped off when external circulation is started.
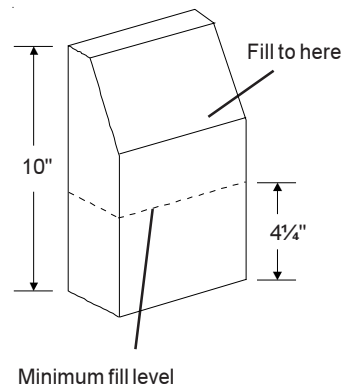
## Reservoir Tank Isolation Valves Open Tank Circulation

Units are equipped with two reservoir tank isolation valves. The valves are used to isolate the tank from the fluid flow. The units are factory shipped with the valves open.

Ensure the valves are open when circulating to a closed system. When operating at low temperatures and circulating to a closed system, should condensation appear on the reservoir, closing the upper valve *after initial start up* will help eliminate it. Leave the valve open for initial start up so process fluid flow through the reservoir can rapidly bleed air from the system. To close the upper valve, open the reservoir access panel and turn the valve 90°.

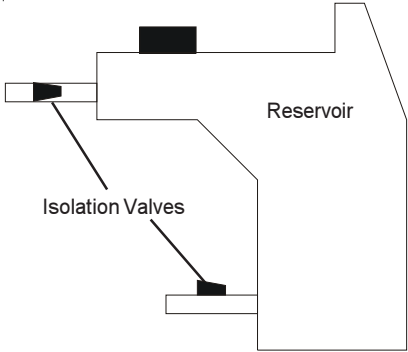Upper Left Side of Unit, Panel Removed


**System temperature changes with the isolation valves closed, and while circulating to a closed system, can create a vacuum and/or high pressure conditions.**

If your application requires the process fluid to flow to an open tank located above the unit, before starting close the valves to prevent fluid siphoning and reservoir overflow. To close the upper valve, open the reservoir access panel, to close the lower valve you must remove the side access panel. Turn the valves 90°.

## Controller

The controller controls temperature using a Proportional-Integral-Derivative (PID) algorithm. It is designed with self-diagnostic features and easy to use operator interface.

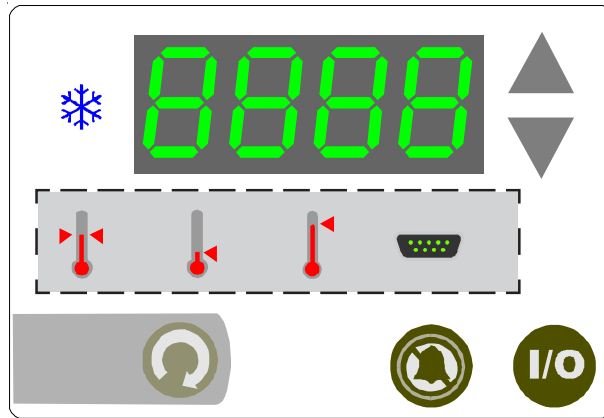Controller

Use this key to toggle the **unit on or off**.

Use this key to **scroll** through the controller's LEDs. It is also used to accept and save changes.

If the alarm sounds, use this key to toggle the **alarm off and on**.

**NOTE:** If the cause of the alarm is cleared but then reoccurs, the alarm will sound again.

Use this key to **increase displayed numerical values**.

Use this key to **decrease displayed numerical values**.

Indicates **refrigeration system status**. It illuminates to indicate the refrigeration system is removing heat from the cooling fluid. As the operating temperature approaches the temperature setpoint, the LED will flash. The indicator is off when heat is not being removed.

**NOTE:** The following indications are not visible until  is depressed.

 Indicates the controller is displaying the **setpoint**. Adjust the display value using the arrow keys.

**NOTE:** You cannot adjust the setpoint closer than 2°C to either of the temperature limits discussed below.


 Indicates the controller is displaying the **low temperature alarm** setting. The indicator flashes and the alarm sounds if this limit is exceeded. Adjust the display value using the arrow keys.

 Indicates the controller is displaying the **high temperature alarm** setting. The indicator flashes and the alarm sounds if this limit is exceeded. Adjust the display value using the arrow keys.

 **(Optional)** Indicates the unit is in the **remote/serial communication** mode of operation. The indicator flashes if the unit is in the remote start/stop mode of operation. Toggle between the local and remote modes using the arrow keys.

**NOTE:** For all four indicators, if neither arrow key is pressed within 10 seconds the display will return to the recirculating temperature display. Save any change using . After pressing an arrow key, if  is not pressed within 60 seconds the display will return to the recirculating temperature display and ignore any change.

## Start Up/Shut Down

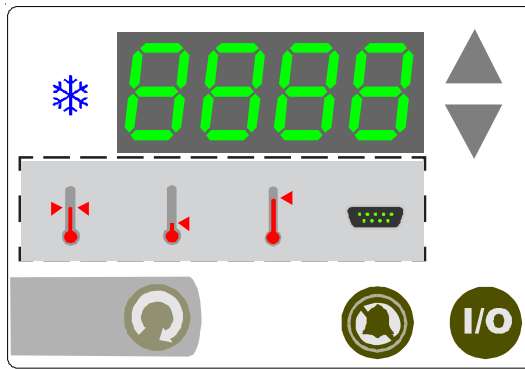Controller

**Pre-Start** Before starting the unit, double check all electrical and plumbing connections. Have extra recirculating fluid on hand. If the unit will not start refer to Section VI Troubleshooting.

**Circuit Breaker** Place the circuit breaker located on the rear of the unit to the up position, the controller will flash and the alarm will momentarily sound.

**The unit has automatic restart. If the unit was shutdown as a result of a power failure and power is restored, it will restart. If it was shut down using the circuit breaker, placing the circuit breaker up will restart the unit.**

**Starting** Press . The controller will do a self-test (the controller will quickly sequence through its LEDs and momentarily sound the alarm) and then display the recirculating fluid temperature. The refrigeration system and the recirculation pump will then start. The RECIRCULATING PRESSURE gauge will display the pump operating pressure. If the pressure needs adjusting, refer to page 40.

**NOTE:** If on start up the unit's recirculating fluid is outside either temperature limit, the unit will operate but the appropriate indicator will flash until the fluid is within the limit.

**Setpoint** To display/change the setpoint press  until  illuminates. The display will flash the current setpoint value. Use the arrow keys to change the value.

**NOTE:** If the arrow keys are not pressed within 10 seconds the display will return to the current reservoir temperature. The controller will not allow you to enter a setpoint closer than 2°C of either temperature alarm setting discussed on the next page. Trying to use a setpoint within 2°C causes the appropriate indicator to flash and sounds the audible alarm twice.

Once the desired setpoint is displayed, press  again to confirm the change. The display will rapidly flash the new value for a short period and then return to the recirculating fluid temperature.

**NOTE:** If the new value is not confirmed within 60 seconds the display will return to the recirculating fluid temperature and ignore any changes.

## Temp Alarms

To display/change the temperature alarm setting press  until either  or  illuminates. The display will flash the current limit value. Use the arrow keys to change the value. The low-end range is 0°C (-15°C for low-temp units) to 30°C. The high-end range is 10°C to 40°C.

**NOTE:** If the arrow keys are not pressed within 10 seconds the display will return to the current reservoir temperature. You cannot set either alarm closer than 2°C of the setpoint.

Once the desired setting is displayed, press  again to confirm the change. The display will rapidly flash the new value for a short period and then return to the recirculating fluid temperature. See Error Codes on page 37.

**NOTE:** If the new value is not confirmed within 60 seconds the display will return to the recirculating fluid temperature and ignore any changes.

## Remote/Serial Comm

**(Optional)** Once your remote/serial comm cable is connected to the unit, enable or disable remote operation or serial communication by pressing  until  illuminates. The display will flash the communication status. Use the up arrow key to enable, the down arrow to disable. See pages 33 and 34 for additional information.

**NOTE:** If the arrow keys are not pressed within 10 seconds the display will return to the current reservoir temperature.

Once the desired status is displayed, press  again to confirm the change. The display will rapidly flash the new status for a short period and then return to the recirculating fluid temperature.

**NOTE:** If the new status is not confirmed within 60 seconds the display will return to the recirculating fluid temperature and ignore any change.

**NOTE:** Enabling serial communications disables the remote start/stop feature. In either mode, the setpoint and alarm values can not be changed using the controller but they can still be viewed.


## External Sensor

**(Optional)** Once the external sensor cable is connected to the unit, enable or disable the sensor by pressing  a fifth time. The display will alternate with rPrb and either  or . If desired, use the arrow keys to change the status and then press  to return to the temperature display. See page 35 for additional information.

## Stopping

To turn the unit off, press .
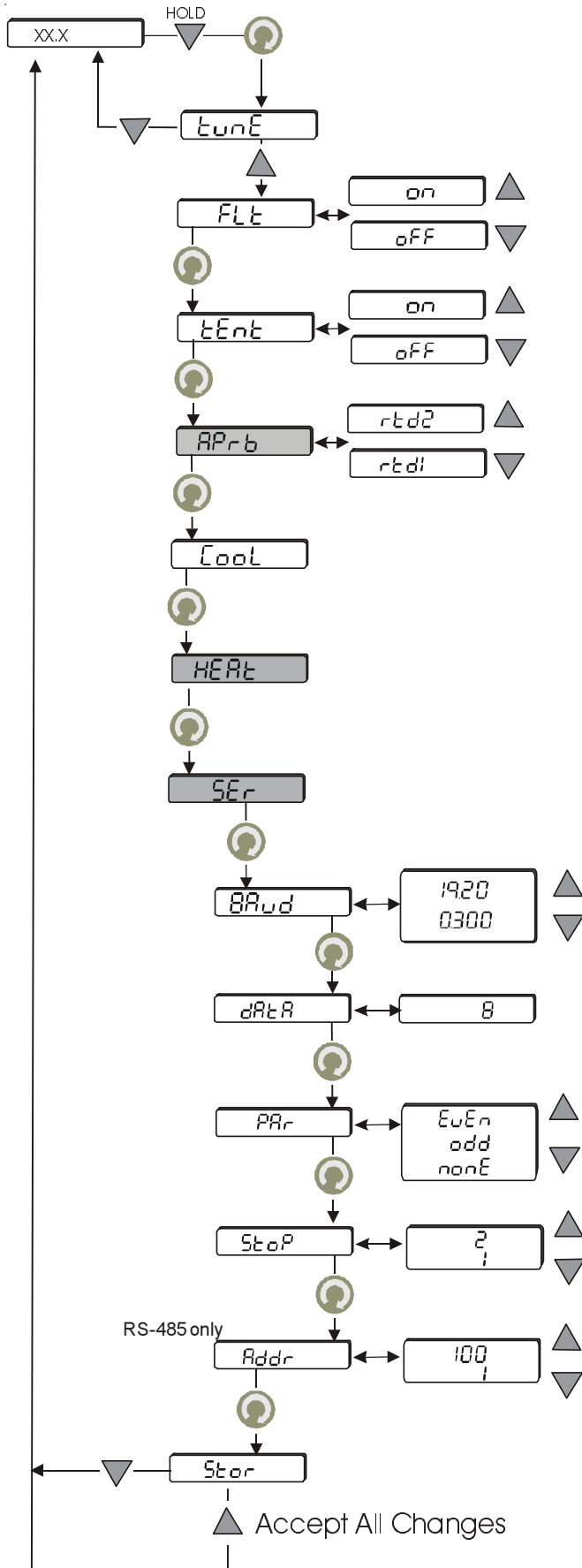
## Setup/Tuning Loop

The Setup/Tuning Loop is used to configure the controller temperature display and various operating parameters. To enter the loop you must be displaying the reservoir fluid temperature and then, while pressing and holding , press . The display will indicate `tunE`. Use the arrow keys to enter/bypass the loop. Once in the loop press  to sequence through it. Use the arrow keys to change any display.

The loop is used to determine how the unit will react when a fault occurs (`FLt`) — either shut down (`on`) or continue to run (`off`). The unit is shipped configured to run.

It is also used to configure the temperature displays to indicate to a tenth of a degree (`tEnt`).

The loop is used to select which temperature sensor the high and low temperature monitors will use (`APrb`). `rtd1` is the internal sensor, `rtd2` is the optional external sensor. `rtd2` is displayed only if the optional external sensor is enabled.

The loop can also be used to set the controller's PID parameters (`Cool HEAt`). The Thermo factory preset values are on page 38.

**Thermo recommends that only a qualified technician change PID parameters. Incorrect PID values will hamper unit performance.**

The heat PID is only applicable to units configured with a heater.

The loop is also used to adjust/verify serial communication parameters (`SEr`). This section of the loop is applicable only to units configured for serial communications capability.

Changes made in the Tuning/Setup Loop take affect when  is pressed at the `StoP` prompt. Pressing  at the prompt aborts *all* changes.

Blank page.

## Section IV Options

### Communication Connector

A female 9-pin D-connector is located on the rear of the unit. The connector is used for either RS-232 or RS-485 serial communication (see Appendix B or C).

When the unit is configured for serial communication the  LED on the controller is illuminated. If the unit is shut off while still configured for serial communication, the temperature display will go blank, but the  LED will still be illuminated.

When the unit is configured for serial communication the unit can be stopped using the local controller. To restart, send another start command.

**NOTE:** If you need to start the unit using the local controller, simultaneously depress and hold both arrow keys for approximately 10 seconds. The display will show the temperature, and the alarm will sound. Use  and scroll to the  LED and turn it OFF. The unit can now be started and operated from the keypad.

Options

#### RS-232 COMM


Pin #	Function
1	No connection
2	TX
3	RX
4	No connection
5	GND = Signal ground
6	No connection
7	No connection
8	No connection
9	No connection

#### RS-485 COMM

Pin #	Function
1	No connection
2	No connection
3	No connection
4	No connection
5	No connection
6	No connection
7	No connection
8	TX
9	RX

TX = Transmitted data from controller  
 RX = Received data to controller.

**Hardware Mating Connector**  
**AMP Part# 745492-2 or equivalent**
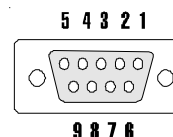
## Remote Start/Stop

**24V**

A male 9-pin D-connector is located on the rear of the unit. The connector is used for either 12 - 24VAC or 12 - 24VDC remote start/stop; and as a status relay dry contact.

To start the unit from a remote location apply the correct voltage to pins 6 & 8. When the unit is started by a remote start/stop voltage  on the controller will flash.

Once started remotely, stop the unit by removing the applied voltage from pins 6 & 8 or by pressing  on the local controller. To restart remotely after stopping locally, remove and reapply the remote start voltage to J101.


### NOTES:

Enabling serial communications disables the remote start/stop feature.

Should the unit experience a high/low temperature warning/fault, a reservoir level warning/fault, or a low pump flow warning/fault, the status relay will de-energize.

### COMM Pin Connection

Pin #	Function
1	No Connection
2	Status, NO
3	Status, NC
4	RTD2A
5	Status, Common

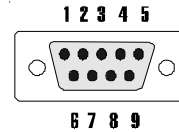**Hardware Mating Connector**  
**AMP Part# 745491-2 or equivalent**

### Remote Start Input:

12 to 24 Volts DC  $\pm 10\%$  maximum input current, 9ma.

12 to 24 Volts AC  $\pm 10\%$  maximum input current, 5ma.

### Absolute maximum input voltages:

40 volts rms. continuous.

56 volts peak.

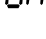
Rating for both normally open and normally closed status relay contacts are:

Maximum rated load - 1.0 Amp, 24 Volts AC/DC, resistive.


Minimum permissible load - 1 ma @ 5 VDC.

## External Temperature Sensor

Connect the external sensor to the male 9-pin D-connector located on the rear of the unit.

The external temperature sensor is enabled using the controller. With the controller displaying the temperature press  five times. The display will alternate with rPrb and either  or . If desired, use the arrow keys to change the status and then press  to return to the temperature display.

If the external sensor is selected the controller will display the temperature at the external sensor. Also, the controller indicates the external sensor is selected by flashing the left-most temperature decimal point.

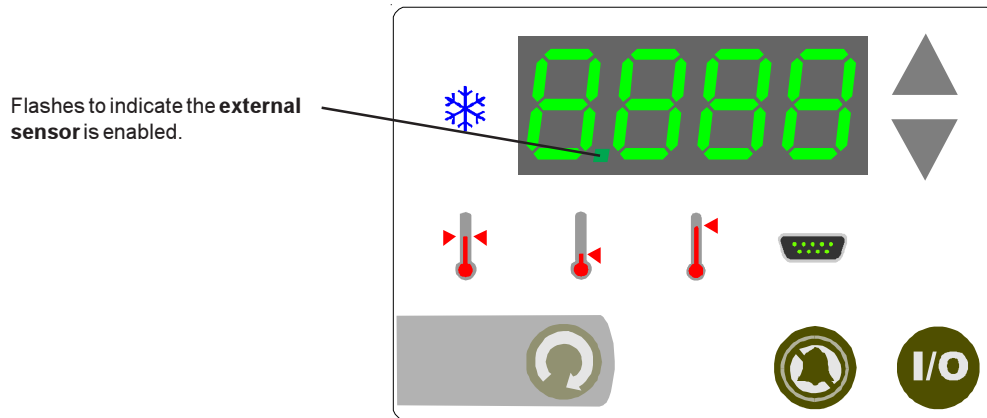The controller continues to use the existing setpoint. When using the external sensor the controller will not allow internal unit temperatures outside 3°C to 37°C for standard temperature range units, -17°C to 37°C for low-temperature range units.

The selected sensor is retained by the controller if it is turned off, or if power is removed.

If there is an external sensor error the unit will automatically return to the internal sensor. An error code will flash and if the fault mode is on the unit will shut down. See Error Codes on page 37. Press  to restart the unit.

The external probe can be used to control the high and low temperature monitors using the Setup/Tuning Loop, see page 31.

**NOTE:** If you need to monitor status while using the external sensor, an optional Y-connector cable is available from Thermo.

## Reservoir Level Warning/Fault

Should the reservoir fluid level drop below normal the mutable alarm will sound, the status relay will de-energize (if the unit is equipped with the communications package), an error message will be displayed, and, depending on the unit's configuration, the unit will continue to run or shut down. The unit is configured to run or shut down using the controller's Setup/Tuning Loop on page 31.

Refill the reservoir and, if necessary, restart the unit.

**NOTE:** The unit will not start without any fluid in the reservoir regardless of how it is configured.

## Low Pump Flow Warning/Fault

Should the pump flow drop below 0.25 gpm (1.0 gpm for units with a CP-55 pump) the mutable alarm will sound, the status relay will de-energize (if the unit is equipped with the communications package), an error message will be displayed, and, depending on the unit's configuration, the unit will continue to run or shut down.

Once the cause of the low flow is identified and corrected restart the unit, if necessary.

The unit is configured to run or shut down using the controller's Setup/Tuning Loop, see page 31.

## Additional Accessories

### External Filtration Package

A partial flow liquid filter available in 5, 10, 25, and 40 micron filter sizes.

### External DEI Package

Partial flow DEI package will maintain between 1 and 3 megohm/cm<sup>2</sup> water resistivity. A cartridge indicator lets you know when the filter needs changing.

### Condenser Air Filter

An easy to install air filter protects the condenser in dusty environments.

### External Strainer Accessory Kit

An easy to install kit to make routine strainer preventative maintenance easier.

## Section V Maintenance

### Error Codes

- Er00* ROM checksum. Clear with  key.
- Er01* Test failure. Locks up the program.
- Er02* Display board failure. Clears when display board is fixed.
- Er03* Critical checksum failure. Clear with  key.
- Er14* Synchronous comm
- Er15* Asynchronous comm
- Er1b* Bad calibration data

### Unusual Hardware Conditions

These errors will flash on the display and cannot be cleared. These are internal controller problems.

- Er04* through *Er13* - Interrupt errors during runtime.
- ConF* BOM invalid

### Functional/ Machine errors

These errors will clear themselves once the problem disappears.

Display	Indication
<i>LoT</i>	Low temp setpoint warning/fault
<i>AdD</i>	Low level warning, fluid below minimum fill level, see page 25
<i>HiT</i>	High temp setpoint warning/fault
<i>Er23</i>	Shorted external temperature sensor (rdt2)
<i>Er24</i>	Open external temperature sensor (rdt2)
<i>Er25</i>	Shorted internal temperature sensor (rdt1)
<i>Er2b</i>	Open internal temperature sensor (rdt1)
<i>Er33</i>	Reservoir fluid below 1°C, standard-temp units only
<i>LoFL</i>	Low process flow warning/fault
<i>LLF</i>	Fluid low level fault - unit shuts down

Unit reaction to warning/fault errors depend on how the unit is configured, see Setup/Tuning Loop on page 31. The unit is shipped configured to continue running. If any other code appears contact our customer service.

## Displaying Software Version Number

The controller can display the installed software version number. For example, for a unit with software version **026950.1A**:

1. Unit is running normally and displaying recirculating fluid temperature.
2. Press and hold  for at least 10 seconds. The display will show the first two digits, for example: *02*.
3. Press .  
The display will show the remaining digits to the left of the decimal, for example: *6950*.
4. Press .  
The display will show the decimal point and the digit to the right of the decimal point, for example, *.1*.
5. Press .  
The display will show the revision letter, as its equivalent number, for example, *1 = A*.
6. Press .  
Disregard this display.
7. Press .  
The display returns to the recirculating fluid temperature.

## Controller PID Values

The factory preset PID values are:

COOL	P	I	D
STD TEMP UNITS	20	0.5	0.0
LOW TEMP UNITS	12	0.5	0.0
HEAT			
ALL UNITS	5	0.5	0.0

**Thermo recommends that only a qualified technician change PID parameters. Incorrect PID values will hamper unit performance.**

## Cleaning

### Reservoir

Periodically inspect the fluid inside the reservoir. If cleaning is necessary, flush the reservoir with a cleaning fluid compatible with the circulating system and the cooling fluid.

The cooling fluid should be replaced periodically. Replacement frequency depends on the operating environment and amount of usage.

**Before changing the cooling fluid ensure it is at a safe handling temperature.**

**NOTE:** It will be necessary to use a wet-vac or siphon to completely drain the reservoir. Use care so as not to damage the reservoir float switch.

See page 25, Filling Requirements, for instructions on cooling fluid replacement.

### Condenser

For proper operation, the unit needs to pull substantial amounts of air through a condenser. A build up of dust or debris on the fins of the condenser will lead to a loss of cooling capacity. Optional air filters are available, contact our Service Department. See Preface, After-Sale Support.

The lower front of the unit has a one-piece grille assembly. Using your hands gently pry the assembly off. Use care not to scratch the paint.

Periodic vacuuming of the condenser fins is necessary. The cleaning frequency depends on the operating environment. After initial installation we recommend a monthly visual inspection of the condenser. After several months, the cleaning frequency will be established.

**Use care cleaning the condenser fins, they can be easily bent.**

### Algae

To restrict the growth of algae in the reservoir, we recommend keeping the reservoir cover in place and using opaque circulation lines. This will eliminate the entrance of light required for the growth of most common algae.

We recommend the use of Chloramine-T, 1 gram per 3.5 liters. Other algaecides can be harmful to the unit's internal components. Contact Thermo for additional information.

### Leaks

A very small leak between the pump and motor is normal. If this is a new pump, allow it to run for 24 hours and reinspect. A continuous leak will require pump replacement.

## Phase Rotation Interlock

Three phase units have a rotation interlock which prevents the unit from starting if the phase rotation is wrong. If the unit will not start, see Section VII, Troubleshooting. If the options in the checklist are not applicable, the problem may be phase rotation.

The interlock is located in the power box. Remove the unit's top and left-side panel to access it. The interlock has a green light which illuminates when the phase is correct.

If the phase is incorrect, disconnect the unit from its power source. Reverse any two line conductors on the line side of the main circuit breaker.

**Never remove the green ground wire.**

Reconnect the unit to its power source. If the unit will not start, contact our Service Department.

## PD Pressure Relief Valve

Units with PD pumps have a Pressure Relief Valve. Refer to the name plate label on the rear of the unit and the BOM decoder on page 40 to identify the type of pump. The valve is used to adjust the unit's fluid pressure.

**NOTE:** The valve is factory preset for the most common applications and normally requires no further adjustment. The valve is factory preset not to exceed  $80 \pm 5$  psi ( $550 \pm 1$  kPa).

**NOTE:** PD pumps are capable of 110 psi maximum output pressure. Ensure your plumbing is rated to withstand 110 psi at the highest operating temperatures.

System back pressure is installation and application dependent. Due to internal Merlin back pressure, the minimum pressure setting is 10 psi for a PD-2 pump, and 6 psi for a PD-1 (these settings prohibit external flow from the Merlin).

Before adjusting the valve turn the unit off. Locate the circular relief valve opening on the rear of the unit. Turn the threaded stem fully clockwise.

If the unit is not plumbed to an application, install a loop of hose equipped with a shut-off valve between the supply and return fittings.

Turn the unit on. Ensure that there is back pressure on the system. Use the pressure gauge to read the relief valve setting.

Back out the threaded stem on the relief valve counterclockwise. Continue until the gauge indicates 80 psi (550 kPa) or the desired setting.

**NOTE:** The relief valve may drip if the threaded stem is backed out too far.

## Pump Strainer

Units with PD pumps, and units with CP-55 pumps that are equipped with reservoir level/low pump flow switches, have a strainer. Refer to the name plate label on the rear of the unit and the BOM decoder on page 42 to identify the type of pump and configuration.

If debris is in the system, the strainer will prevent the material from being drawn into the flow switch and the pump, damaging the pump vanes. A clogged strainer will also cause increased discharge pump pressures.

After initial installation, the strainer may become clogged with debris and scale. Therefore, **the strainer must be cleaned after the first week of installation**. After this first cleaning, we recommend a monthly visual inspection. After several months, the cleaning frequency will be established.

**Before cleaning the strainer ensure the cooling fluid is at a safe handling temperature. Disconnect the power cord from the power source.**


We also recommend draining the reservoir.

The strainer is in the pump return line, the exact location depends on the unit's configuration. For most configurations it is located behind the rear panel. For units configured with low level and low temperature switches, it is behind the right-side panel. **NOTE:** A small amount of fluid will come out of the internal strainer housing when the fitting is removed. Have a rag handy to put under the housing before removing the fitting.

Remove the panel from the unit. Use one wrench to support the filter housing, and a second wrench to unscrew the fitting,  $1\frac{3}{8}$  inch for CP pumps, or  $\frac{3}{4}$  inch for PD pumps. Remove the screen. **NOTE:** Some units with  $\frac{3}{4}$  inch fittings also have a brass end cap, see illustration below.

Clean the screen by rinsing it with water. When the screen is clean, replace it and the cap. Tighten the fitting and replace the panel. Refer to Section III for instructions on replacing the cooling fluid.


**NOTE:** An external strainer accessory kit is available to make routine strainer preventative maintenance easier. The part number is 62700000010.
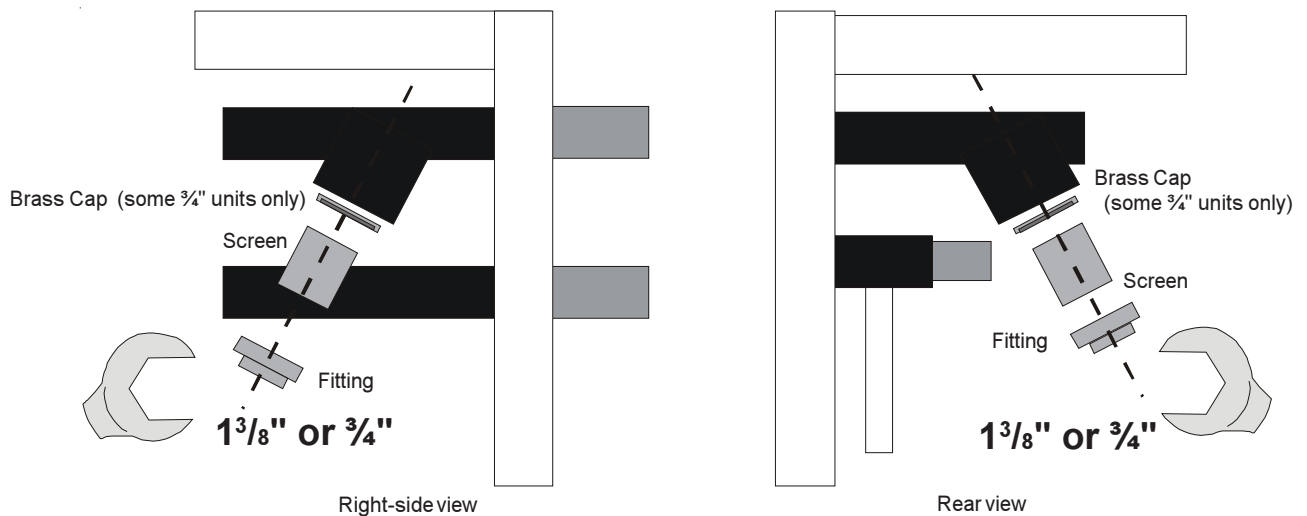
## PD Pump Motor Lubrication

Some PD pumps use sleeve type bearings with large reservoirs. Oiling instructions are generally posted on each motor. In the absence of instructions, add approximately 30 to 35 drops of SAE 20 (142 CS viscosity) non-detergent oil (i.e. 3-IN-ONE® Tool Oil) to each bearing on the following schedule:

Duty Cycle	Oiling Frequency
Continuous	Once every year
Intermittent	Once every 2 years
Occasional	Once every 5 years

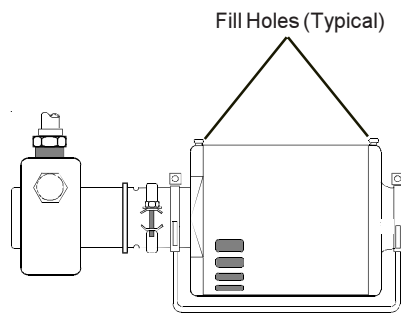
### Storage

If the unit is to be transported and/or stored in cold temperatures it needs to be drained and then flushed with a 50/50 glycol/water mixture.

### BOM Decoder

The Bill of Material (BOM) number helps identify the configuration of your unit. The number is printed on the label located on the rear of the unit.

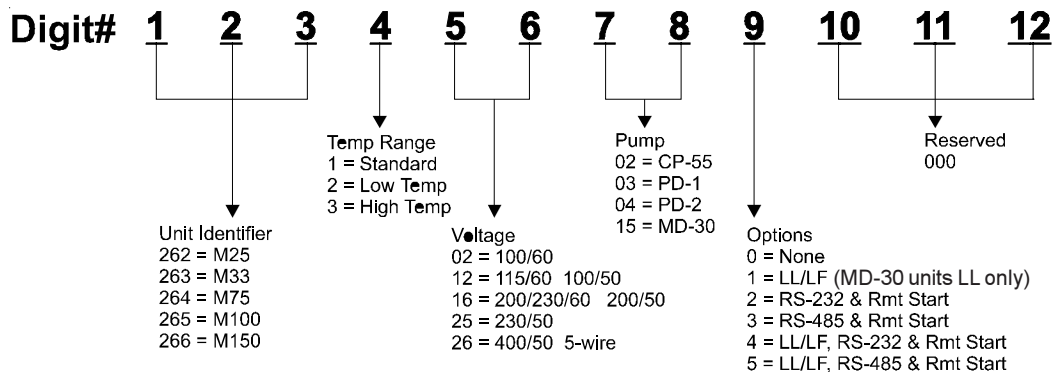
# Temperature Sensor Calibration

## Internal (rtd1)

If the temperature fluid temperature display disagrees with your reference thermometer, the internal temperature sensor (rtd1) may need calibration.

Do not pick calibration points that are outside the safe operating limits of the fluid in your application. For example with water, 35°C and 5°C would be typical calibration points.

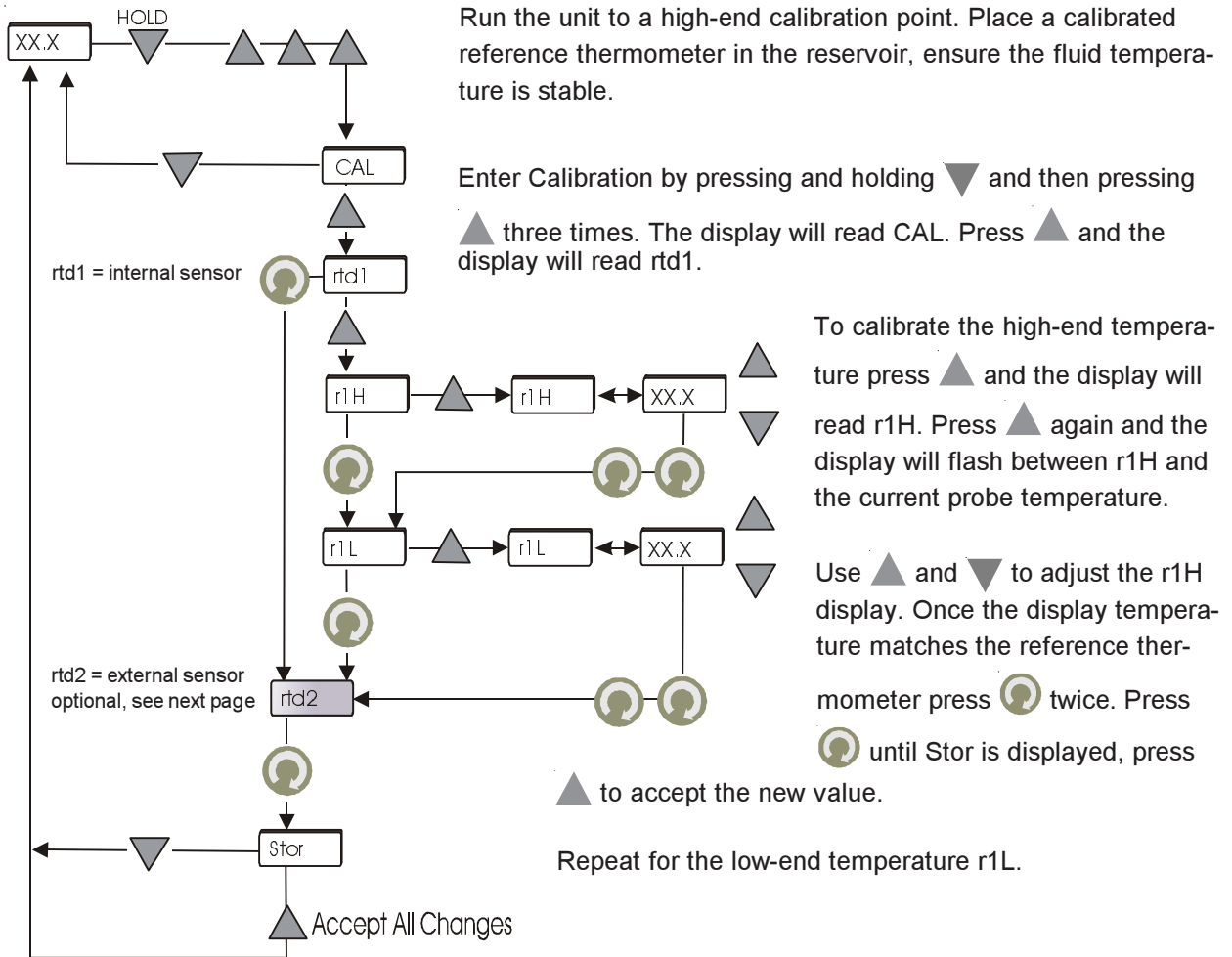
**NOTE:** Both the high and low temperatures must be entered for a valid calibration. If power is lost before the procedure is complete, critical data needed for the calculation of calibration parameters will be lost. After pressing ▲ at the Stor prompt wait several seconds before proceeding to ensure that a bad calibration message does not appear. Premature use of the keypad after pressing ▲ may cancel the bad calibration error message.

### External (rtd2)

When using any external sensor for the first time, Thermo recommends a two-point calibration. This calibration will only affect the temperature read by the external sensor.

Do not pick calibration points that are outside the safe operating limits of the fluid in your application. For example with water, 35°C and 5°C would be typical calibration points.

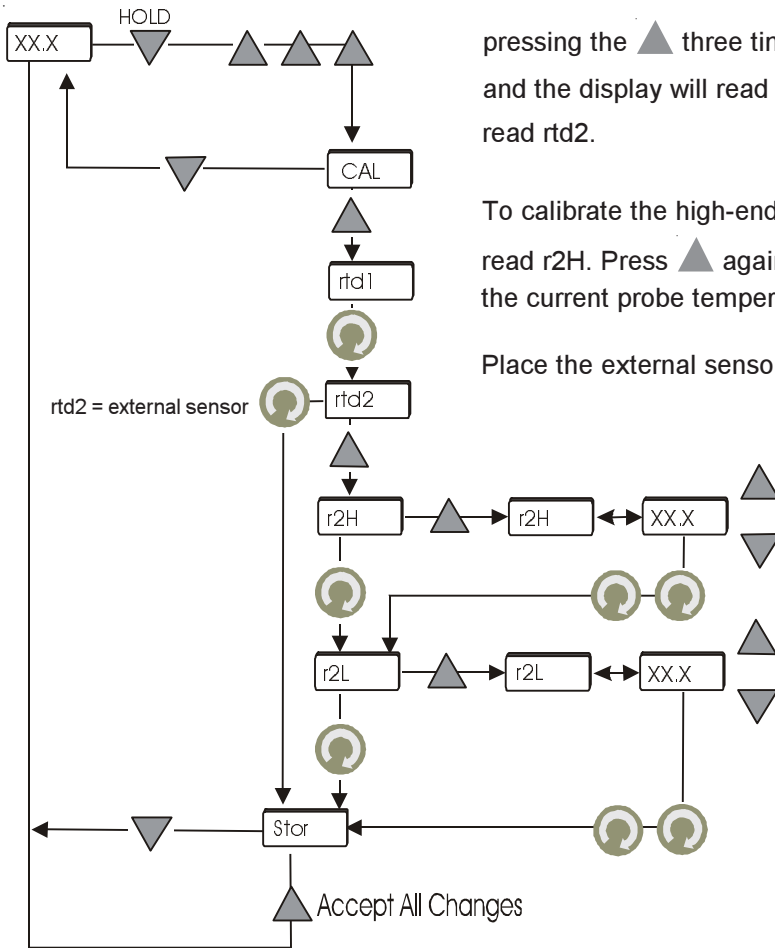Enter the Calibration Loop by pressing and holding ▼ and then pressing the ▲ three times. The display will read CAL. Press ▲ and the display will read rtd1, then press → and the display will read rtd2.

To calibrate the high-end temperature press ▲ and the display will read r2H. Press ▲ again and the display will flash between r2H and the current probe temperature.

Place the external sensor and a calibrated reference thermometer in a source containing fluid that is just below the high-end temperature range of the unit, ensure the fluid temperature is stable. Use ▲ and ▼ to adjust the temperature display to match the reference thermometer.

Once the display temperature matches the reference thermometer press → twice. Press → until Stor is displayed, pressing ▲ accepts the new value.

Repeat for the low-end temperature, r2L. Use a fluid that is just above the low-end temperature range of the unit.

**NOTE:** Both the high and low temperatures must be entered for a valid calibration. If power is lost before the procedure is complete, critical data needed for the calculation of calibration parameters will be lost. After pressing ▲ at the Stor prompt wait several seconds before proceeding to ensure that a bad calibration message does not appear. Premature use of the keypad after pressing ▲ may cancel the bad calibration error message.

## Section VI Troubleshooting

### Checklist

#### Unit will not start/shuts down

Check the line cord; ensure it is plugged in.

Check the position of the circuit breaker on the rear of the unit.

Check the controller for error codes, see page 37.

**NOTE:** On units with a Low Flow Switch and configured to shut down with a low flow fault (L<sub>o</sub>FL), several starting attempts may be necessary.

On units with a Reservoir Level Switch check the reservoir level.

Optional serial communications or remote start/stop enabled.

Make sure power source voltage is within the unit rated voltage  $\pm 10\%$ .

400 Volt, 50 Hertz units are equipped with a phase monitor. Ensure proper phase to the unit. See Phase Rotation Interlock on page 40.

If  is lit, the unit was shut off using serial communications. Start the unit using serial communications. To start the unit using the local controller press and hold both arrow keys for 10 seconds, the controller will power. Turn serial communications off using the procedure on page 33, and then press  to turn off the controller. Press  again to start the unit.

#### Unit will not circulate fluid

Check the reservoir level. Fill, if necessary.

Check the position of the reservoir tank valves. See Reservoir Tank Isolation Valves/Open Tank Circulation.

Check the instrument being cooled for restrictions in the cooling line.

Check the pump strainer (units with PD pumps).

The pump motor may be too hot. Allow it to cool.

Check the pressure gauge, adjust the relief valve as necessary (units with PD pumps).

On units with MD pumps never completely restrict flow to your application. Dead-heading the pump will damage its coupling and require pump replacement.

Continued on next page.

**Inadequate temperature control**

Verify the setpoint.

If the temperature continues to rise, make sure your application's heat load does not exceed the rated specification, see Section II, Specifications.

Make sure the air intake and discharge are not impeded and the ambient temperature does not exceed +35°C.

Make sure the condenser is free of dust and debris, see Section V, Maintenance and Cleaning.

Verify the controller's PID values, see Setup/Tuning Loop.

**Unit runs loudly**

Check the position of the reservoir tank valves. See Reservoir Tank Isolation Valves on page 26..

**Condensation on reservoir**

Check the position of the upper reservoir tank valve. See Reservoir Tank Isolation Valves on page 26.

**No serial communications (optional)**

All units are tested for serial communications before they leave the factory. Check the controller, ensure serial communication is enabled.

Check all communications commands, they *must be exact*. See Appendix B or C.

Check all wiring for proper connections or possible shorts.

Software to verify serial communication is available from Thermo.

**Service Assistance**

If, after following these troubleshooting steps, your unit fails to operate properly, contact our Service Department for assistance (see Preface, After-sale Support). Before calling, *please* obtain the following information:

- *unit BOM number*
- *unit serial number*
- *software version*
- *voltage of power source*

## Appendix A Water Quality Standard and Recommendations

	Permissible (PPM)	Desirable (PPM)
<b>Microbiologicals</b>		
(algae, bacteria, fungi)	0	0
<b>Inorganic Chemicals</b>		
Calcium	<40	<0.6
Chloride	<250	<25
Copper	<1.3	<1.0
Iron	<0.3	<0.1
Lead	<0.015	0
Magnesium	<12	<0.1
Manganese	<0.05	<0.03
Nitrates\Nitrites	<10 as N	0
Potassium	<20	<0.3
Silicate	<25	<1.0
Sodium	<20	<0.3
Sulfate	<250	<50
Hardness	<17	<0.05
Total Dissolved Solids	<50	<10
<b>Other Parameters</b>		
pH	6.5-8.5	7-8
Resistivity	0.01*	0.05-0.1*

\* Megohm-Cm (Compensated to 25°C)


Unfavorably high total ionized solids (TIS) can accelerate the rate of galvanic corrosion. These contaminants can function as electrolytes which increase the potential for galvanic cell corrosion and lead to localized corrosion such as pitting which can be observed at the studs and on the outside surface of cooling coils. Eventually, the pitting will become so extensive that the coil will leak refrigerant into the water reservoir.

For example, raw water in the U.S. averages 171 ppm (of NaCl). The recommended level for use in a water system is between 0.5 to 5.0 ppm (of NaCl).

Recommendation: Initially fill the tank with distilled/deionized water. Do not use untreated tap water as the total ionized solids level may be too high.

Maintain this water quality at a resistivity of between 1 to 10 megohm-cm (compensated to 25°C) by using a purification system. Although the initial fill may be as high as 10 megohm-cm (compensated to 25°C), the desired level for long time usage is 1 to 3 megohm-cm (compensated to 25°C).

The above two recommendations will reduce the electrolytic potential of the water and prevent or reduce the galvanic corrosion observed.
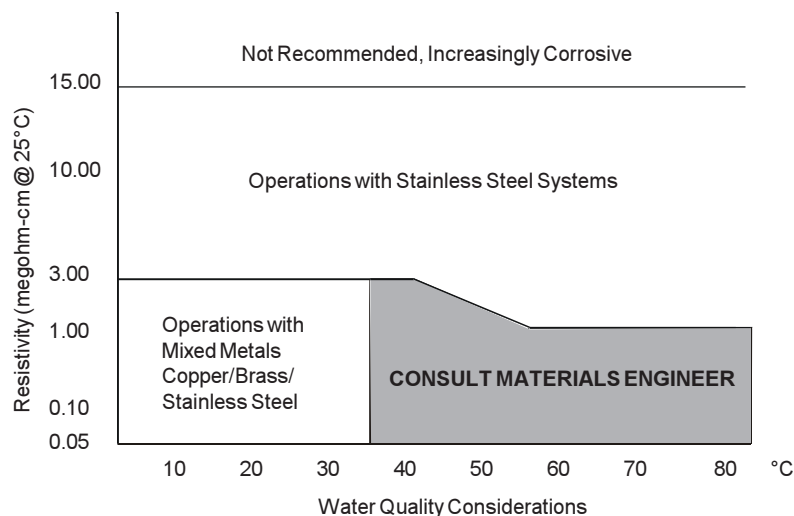
## Merlin Compatibility with Recommended Fluids

### **Filtered/Singled Distilled water, +5°C to +35°C**

This fluid is recommended primarily because it has all microorganisms that cause biological fouling removed through vaporizing and condensing the water. However, distilled water does not remain pure for very long when exposed to the atmosphere. Air-born spores can contaminate the water and activate alga growth. Chloramine-T is a compatible algaecide that can be used to combat growth but a more effective maintenance plan would include switching out the fluid with newly distilled water every six month. The particulates that have been filtered out in the process are also preventative in keeping the system “clean” of contaminants. One thing to note is that distilling water that contains an additive could increase the concentration of that additive in the water.

### **50/50 Uninhibited Ethylene Glycol/Water (-15°C to + 35°C)**

Ethylene glycol is used to depress the freezing point of water as a coolant. We recommend not using the uninhibited (no corrosion additives) ethylene glycol. It is more corrosive to copper than plain water so therefor it should not be recommended unless it is required for the application.

### **50/50 Inhibited Ethylene Glycol/ Water (-15°C to +35°C)**

Inhibited glycol can be used to increase the operating temperature range of the fluid but should not be used as a “pre-mixed anticorrosive” solution. Dow uses a pH standard of 8 to base when the fluid has become corrosive. Dowtherm is an ethylene based product that contains dipotassium phosphates in a 4% concentration. The recommended use of Dowtherm is mixing with distilled or deionized water or water that contains less than 25 ppm chloride and sulfate and less than 100 ppm total hardness of CaCO<sub>3</sub>. The general term, inhibited glycol/water, almost too close to meaning inhibited water. Inhibited water can have many types of additives including chromate that will foul our cooling system very quickly. Some inhibitor additives can release the bonding agent in the carbon graphite in the PD2 pumps and therefor are incompatible, such as Sodium Hydroxide.

### **50/50 Uninhibited Propylene Water (-15°C to +35°C)**

Although the use of this glycol similar to ethylene glycol, propylene glycol is considered “safe” to use in the food industry. Propylene is less dense than ethylene and therefor will have a tendency to weep more through mechanical seals.

### **50/50 Inhibited Propylene/ Water (-15°C to +35°C)**

Same issues as with uninhibited propylene and inhibited ethylene glycol.

### **Deionized water (1-3 megohm, compensated, +5°C to +35°C)**

This water has the ions controlled so that they will not conduct and cause galvanic corrosion between dissimilar metals. Deionized water is aggressive to metal when it is too pure since this is not the normal state of water. The leaching of metallic ions from the metal surface is seen in pitting.

**NEVER use flammable or corrosive fluids with this unit. Do not use automotive anti-freeze. Commercial antifreeze contains silicates that can damage the pump seals. Use of automotive antifreeze will void the manufacturer’s warranty.**

## Wetted Materials

### PD PUMPS (PD1 Part Number 008816, PD2 Part Number 008806):

Carbon Graphite (Graphitar Type B)

Brass

303 Stainless Steel

Ceramic (85% Aluminum Oxide and Balance is Considered Glass)

Fluorocarbon (Viton equivalent)

Ultem (relief Valve Material)

### MD-30 (Part Number 026922):

Polypropylene

Ceramic

Magnet

Buna-N

### CP-55 PUMP (Part Number 026652):

Stainless Steel

Viton

### Evaporator

#### (Part Numbers 026964/026945/026880):

AISI 316 Stainless Steel

Copper Brazing

### Tank Fittings:

Brass

### Strainer (Part Number 009334):

Bronze

304 Stainless Steel Mesh

### Tank Cap (Part Number 026992):

Nylon

### Tank Body (Part Number 026933):

Polyethylene (LLPDE, Mobile NRP=135)

### Hose Assemblies/Plumbing:

Buna-N (Parker Push-Lok O-ring)

Nitrile based synthetic rubber (Parker Push-Lok hoses) Part Number 026931

Brass Part Number 032688

Copper Part Number 003752

### Ball Valve (Part Number 026981):

Viton

Chrome plated brass body and ball

PTFE

### Sealants:

FasSeal

Teflon tape

Teflon paste

Nickel Impregnated Teflon tape

Rector-Seal Tplus2

### Low Level/Low Flow Switches (Optiona1):

Polysulphone (Part Number 026900) Low Level

Brass (Part Number 026901) Low Flow

Blank page.

## Appendix B RS-232 NC Serial Communications Protocol


**NOTE:** This appendix assumes you have a basic understanding of communications protocols.

All data is sent and received in binary form, do not use ASCII. In the following pages the binary data is represented in hexadecimal (hex) format.

The NC Serial Communications Protocol is based on a master-slave model. The master is a host computer, while the slave is the chiller's controller. Only the master can initiate a communications transaction (half-duplex). The slave ends the transaction by responding to the master's query. The protocol uses an RS-232 serial interface with the default parameters: 9600 baud, 8 data bits, 1 stop bit, no parity and no handshaking.

**NOTE:** Before the unit will communicate, RS-232 must be turned on using the controller. Ensure the  indicator is illuminated.

The unit can be controlled through your computer's serial port by using the unit's standard 9-pin RS-232 connection. Data read of the serial port connects to the data transmit (pin 2) of the unit. Data transmit of the serial port connects to data read (pin 3) of the unit.

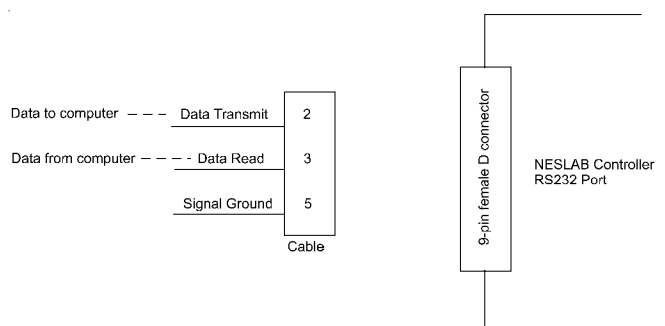Communication cables are available from Thermo. Contact our sales department for additional information.

All commands must be entered in the exact format shown in the tables on the following pages. The tables show all commands available, their format and responses. Controller responses are either the requested data or an error message. The controller response *must* be received before the host sends the next command.

The host sends a command embedded in a single communications packet, then waits for the controller's response. If the command is not understood or the checksums do not agree, the controller responds with an error command. Otherwise, the controller responds with the requested data. If the controller fails to respond within 1 second, the host should resend the command.

# NC Serial Communications Protocol

**NOTE:** All byte values are shown in hex, hex represents the binary values that must be sent to the chiller. **Do not use ASCII.**

The framing of the communications packet in both directions is:

Checksum region								
Lead char CA	Addr-MSB 0	Addr-LSB 1	Command	n d-bytes	d-byte 1	...	d-byte n	Checksum

<i>Lead char</i>	CA (hex).
<i>Addr-msb</i>	Device address is 1
<i>Addr-lsb</i>	Most significant byte of device address is 0.
<i>Command</i>	Least significant byte of device address is 1.
<i>n d-bytes</i>	Command byte (see Table 1).
<i>d-byte 1</i>	Number of data bytes to follow (00 to 03 hex).
...	1 <sup>st</sup> data byte (the qualifier byte is considered a data byte).
<i>d-byte n</i>	...
<i>Checksum</i>	n <sup>th</sup> data byte.
	Bitwise inversion of the 1 byte sum of bytes beginning with the most significant address byte and ending with the byte preceding the checksum. (To perform a bitwise inversion, "exclusive OR" the one byte sum with FF hex.)


The master requests information by sending one of the Read Functions as shown in Table 1. Since no data is sent to the chiller during a read request, the master uses 00 for the number of data bytes following the command byte.

**Appendix B**

The chiller will respond to a Read Function by echoing the lead character, address, and command byte, followed by the requested data and checksum. When the chiller sends data, a qualifier byte is sent first, followed by a two byte signed integer (16 bit, MSB sent first). The qualifier byte indicates the precision and units of measure for the requested data as detailed in Table 2.


As an example, the master requests to read internal temperature by sending:

If the temperature is -12°C, the unit would reply:

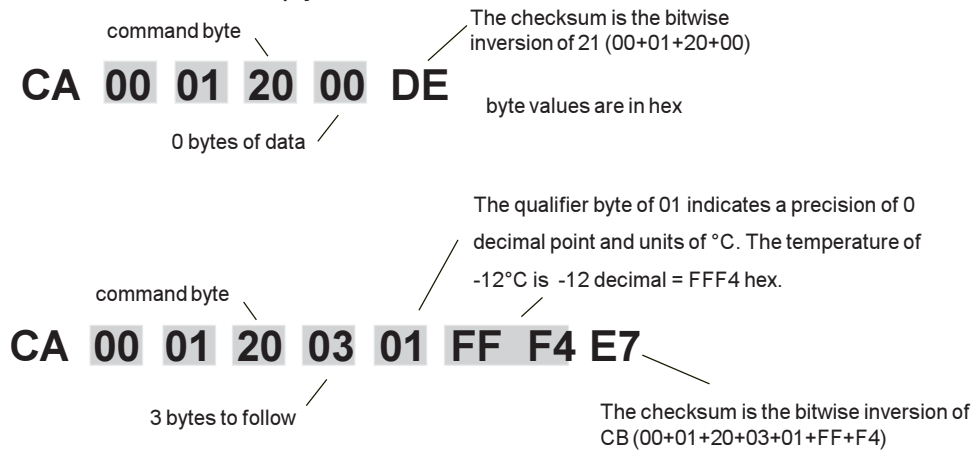
The master sets parameters in the chiller by sending one of the Set Functions as shown in Table 1. The master does not send a qualifier byte in the data field. The master should be preprogrammed to send the correct precision and units (it could also read the parameter of interest first to decode the correct precision and units needed).

For example, if the master wants to set the setpoint to 30°C, it would send :

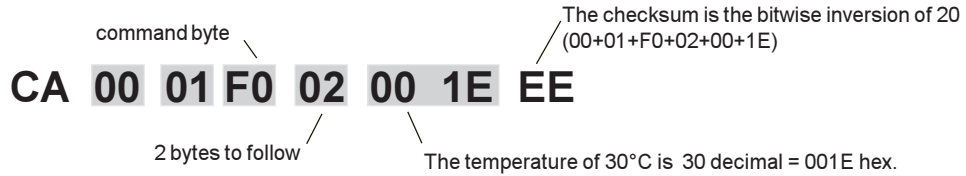The unit responds:

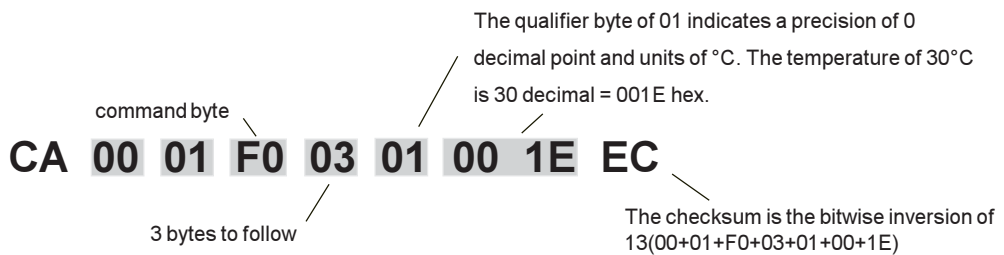**Table 1 Commands (All bytes are in hex)**

FUNCTION	MASTER SENDS	UNIT RESPONDS
Read Acknowledge	CA 00 01 <b>00</b> 00 FE	CA 00 01 <b>00</b> 02 v1 v2 cs
Read Status*	CA 00 01 <b>09</b> 00 F5	CA 00 01 <b>09</b> 02 d1 d2 cs
Error		CA 00 01 <b>0F</b> 02 en ed cs
Read Internal Temperature (RTD1)	CA 00 01 <b>20</b> 00 DE	CA 00 01 <b>20</b> 03 qb d1 d2 cs
Read External Temperature(RTD2)	CA 00 01 <b>21</b> 00 DD	CA 00 01 <b>21</b> 03 qb d1 d2 cs
Read Setpoint (control point)	CA 00 01 <b>70</b> 00 8E	CA 00 01 <b>70</b> 03 qb d1 d2 cs
Read Low Temperature Limit	CA 00 01 <b>40</b> 00 BE	CA 00 01 <b>40</b> 03 qb d1 d2 cs
Read High Temperature Limit	CA 00 01 <b>60</b> 00 9E	CA 00 01 <b>60</b> 03 qb d1 d2 cs
Read Cool Proportional Band (P)	CA 00 01 <b>74</b> 00 8A	CA 00 01 <b>74</b> 03 qb d1 d2 cs
Read Cool Integral (I)	CA 00 01 <b>75</b> 00 89	CA 00 01 <b>75</b> 03 qb d1 d2 cs
Read Cool Derivative (D)	CA 00 01 <b>76</b> 00 88	CA 00 01 <b>76</b> 03 qb d1 d2 cs
Read Heat Proportional Band (P)	CA 00 01 <b>71</b> 00 8D	CA 00 01 <b>71</b> 03 qb d1 d2 cs
Read Heat Integral (I)	CA 00 01 <b>72</b> 00 8C	CA 00 01 <b>72</b> 03 qb d1 d2 cs
Read Heat Derivative (D)	CA 00 01 <b>73</b> 00 8B	CA 00 01 <b>73</b> 03 qb d1 d2 cs

command bytes shown in **bold**  
 qb = qualifier byte, see Table 2  
 d1,d2 = 16 bit signed integer of the value being  
 checksent or received  
 cs = the checksum of the string (see text)  
 \* See Read Status, see Table 3

v = protocol version  
 ed = echoed command byte  
 en 01=Bad command, 02=Bad  
 xx = no valid data, include in checksum

**Table 1 Commands (continued) (All bytes are in hex)**

FUNCTION	MASTER SENDS	UNIT RESPONDS
Set Setpoint (control point)*	CA 00 01 <b>F0</b> 02 d1 d2 cs	CA 00 01 <b>F0</b> 03 qb d1 d2 cs
Set Low Temp	CA 00 01 <b>C0</b> 02 d1 d2 cs	CA 00 01 <b>C0</b> 03 qb d1 d2 cs
Set High Temp	CA 00 01 <b>E0</b> 02 d1 d2 cs	CA 00 01 <b>E0</b> 03 qb d1 d2 cs
Set Cool Proportional Band (P = 0.1 - 99.9)	CA 00 01 <b>F4</b> 02 d1 d2 cs	CA 00 01 <b>F4</b> 03 qb d1 d2 cs
Set Cool Integral (I = 0 - 9.99)	CA 00 01 <b>F5</b> 02 d1 d2 cs	CA 00 01 <b>F5</b> 03 qb d1 d2 cs
Set Cool Derivative (D = 0 - 5.0)	CA 00 01 <b>F6</b> 02 d1 d2 cs	CA 00 01 <b>F6</b> 03 qb d1 d2 cs
Set Heat Proportional Band (P = 0.1 - 99.9)	CA 00 01 <b>F1</b> 02 d1 d2 cs	CA 00 01 <b>F1</b> 03 qb d1 d2 cs
Set Heat Integral (I = 0 - 9.99)	CA 00 01 <b>F2</b> 02 d1 d2 cs	CA 00 01 <b>F2</b> 03 qb d1 d2 cs
Set Heat Derivative (D = 0 - 5.0)	CA 00 01 <b>F3</b> 02 d1 d2 cs	CA 00 01 <b>F3</b> 03 qb d1 d2 cs

**Set On/Off Array**

The Set On/Off Array command is used to set and/or request unit settings. Depending on the software version, the array consists of up to 5 data bytes, d1 - d5. Each byte represents the state of one setting, see \*\* below. The master can send the array with a value of 0 (turn off/disable), 1 (turn on/enable), or 2 (do not change). 0 or 1 causes the unit to change the setting, 2 only requests the state of a setting. The slave returns values of 0 (off/disabled) or 1 (on/enabled).

(SW Version 02950.1A-E)	CA 00 01 <b>81</b> 01 d1 cs	CA 00 01 <b>81</b> 01 d1 cs
(SW Version 02950.1F-H)	CA 00 01 <b>81</b> 04 d1... d4 cs	CA 00 01 <b>81</b> 04 d1...d4 cs
(SW Version 02950.1J-K)	CA 00 01 <b>81</b> 05 d1... d5 cs	CA 00 01 <b>81</b> 05 d1...d5 cs
(SW Version 02950.1L- )	Any of the above formats	
Sample - turn unit on (02950.1J-K)	CA 00 01 <b>81</b> 05 01 02 02 02 02 6F	

\* = limited to the range of the chiller

command bytes shown in **bold**

qb = qualifier byte, see Table 2

d1, d2 = 16 bit signed integer of the value being sent or received

cs = the checksum of the string (see text)

\*\* d1 - unit off = 0, on =1

d2 - external sensor disabled = 0, enabled = 1

d3 - fault mode disabled = 0, enabled = 1

d4 - tenths display disabled = 0, enabled = 1

d5 - Alarms on internal = 0, or external sensor =1, (1 is invalid when d2=0, external sensor disabled)

**Table 2**

QUALIFIER BYTE	QUALIFIER BYTE
00 0 precision, no units of measure	01 precision, °C
10 1 precision, no units of measure	11 precision, °C
20 2 precision, no units of measure	

Example: The temperature of 45.6°C would be represented by the qualifier 11 hex, followed by the 2 bytes 01 C8 hex (456 decimal).

**Table 3**

**READ STATUS**

BIT	d1	d2 <sup>t</sup>
b.7 = 1	External Temp Sensor Enabled	Reserved
b.6 = 1	Level 1 Warning/Auto Refill	Reserved
b.5 = 1	Low Flow Warning	Internal Temp Sensor Fault
b.4 = 1	Low Level 2 Warning	External Temp Sensor Fault
b.3 = 1	High or Low Temp Warning	High Temp/Fixed High Temp Fault
b.2 = 1	High or Low Temp Bypass	Low Temp/Fixed Low Temp Fault
b.1 = 1	Unit Faulted*	Low Flow Fault
b.0 = 1	Unit Running*	Low Level 2 Fault

<sup>t</sup> d2 is the detail of fault d1-b.1

\* mutually exclusive

## Appendix C NC RS-485 Serial Communications Protocol

**NOTE:** This appendix assumes you have a basic understanding of communications protocols.

All data is sent and received in binary form, do not use ASCII. In the following pages the binary data is represented in hexadecimal (hex) format.

The NC Serial Communications Protocol is based on a master-slave model. The master is a host computer, while the slave is the bath's controller. Only the master can initiate a communications transaction (half-duplex). The slave ends the transaction by responding to the master's query. The protocol uses an RS-485 serial interface with the default parameters: 9600 baud, 1 start bit, 8 data bits, 1 stop bit, no parity and a selectable address from 1 to 100.

**NOTE:** Before the unit will communicate, RS-485 must be turned on using the controller. Ensure the  indicator is illuminated.

The unit can be controlled through your computer's RS-485 serial port by using the unit's 9-pin RS-485 connection. Communication cables are available from Thermo. Contact our sales department for additional information.

All commands must be entered in the exact format shown in the tables on the following pages. The Tables on the last page of this Appendix show all commands available, their format and responses. Controller responses are either the requested data or an error message. The controller response *must* be received before the host sends the next command.


The host sends a command embedded in a single communications packet, then waits for the controller's response. If the command is not understood or the checksums do not agree, the controller responds with an error command. Otherwise, the controller responds with the requested data. If the controller fails to respond within 1 second, the host should resend the command.

# NC Serial Communications Protocol

**NOTE:** All byte values are shown in hex, hex represents the binary values that must be sent to the bath. **Do not use ASCII.**

The framing of the communications packet in both directions is:

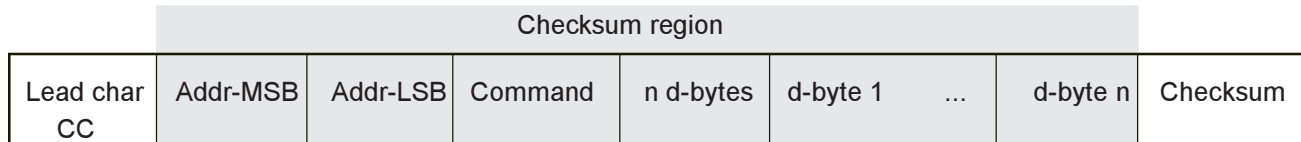
<i>Lead char</i>	CC (hex).
<i>Addr-msb</i>	Most significant byte of device address is 00 hex.
<i>Addr-lsb</i>	Least significant byte of device address in hex (01 to 63 hex).
<i>Command</i>	Command byte (see Table 1).
<i>n d-bytes</i>	Number of data bytes to follow (00 to 03 hex).
<i>d-byte 1</i>	1 <sup>st</sup> data byte (the qualifier byte is considered a data byte).
...	...
<i>d-byte n</i>	n <sup>th</sup> data byte.
<i>Checksum</i>	Bitwise inversion of the 1 byte sum of bytes beginning with the most significant address byte and ending with the byte preceding the checksum. (To perform a bitwise inversion, "exclusive OR" the one byte sum with FF hex.)

The master requests information by sending one of the Read Functions as shown in Table 1. Since no data is sent to the bath during a read request, the master uses 00 for the number of data bytes following the command byte.


The bath will respond to a Read Function by echoing the lead character, address, and command byte, followed by the requested data and checksum. When the bath sends data, a qualifier byte is sent first, followed by a two byte signed integer (16 bit, MSB sent first). The qualifier byte indicates the precision and units of measure for the requested data as detailed in Table 2.

**Appendix C**

As an example, the master requests to read internal temperature unit address 01 by sending:

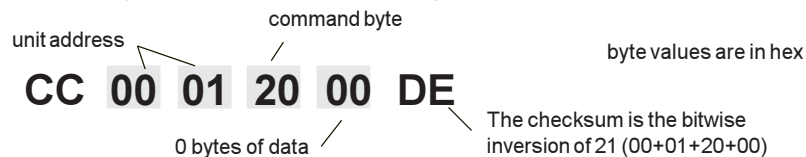
If the temperature is -10.5°C, unit #1 would reply:

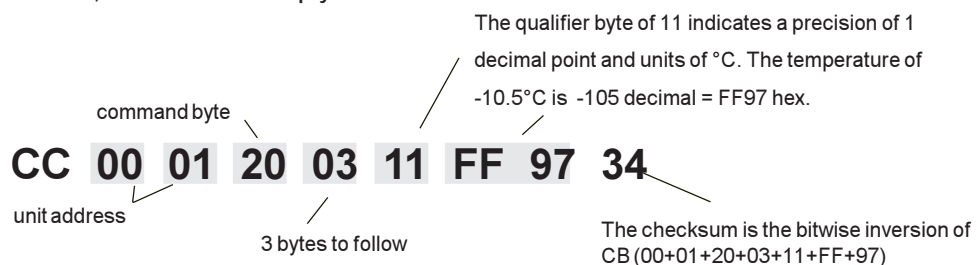
The master sets parameters in the bath by sending one of the Set Functions as shown in Table 1. The master does not send a qualifier byte in the data field. The master should be preprogrammed to send the correct precision and units (it could also read the parameter of interest first to decode the correct precision and units needed).

For example, if the master wants to set the setpoint to 30°C for unit #3, it would send :

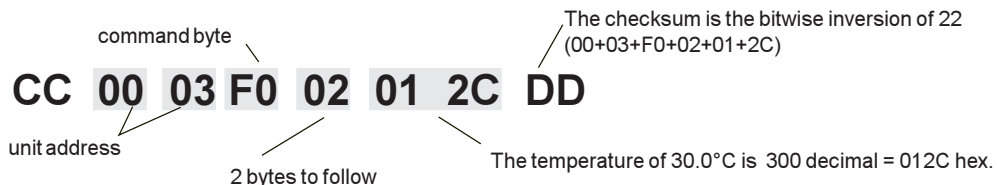Unit #3 responds:

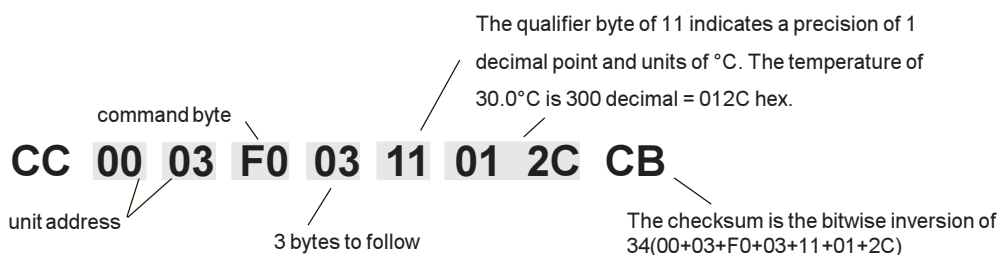**Table 1 Commands (All bytes are in hex)**

FUNCTION	MASTER SENDS	UNIT RESPONDS
<b>READ</b>		
Read Acknowledge	CC 00 01 <b>00</b> 00 FE	CC 00 01 <b>00</b> 02 v1 v2 cs
Read Status*	CC 00 01 <b>09</b> 00 F5	CC 00 01 <b>09</b> 02 d1 d2 cs
Error		CC 00 01 <b>0F</b> 02 en ed cs
Read Internal Temp (RTD1)	CC 00 01 <b>20</b> 00 DE	CC 00 01 <b>20</b> 03 qb d1 d2 cs
Read External Temp (RTD2)	CC 00 01 <b>21</b> 00 DD	CC 00 01 <b>21</b> 03 qb d1 d2 cs
Read External Temp (RTD2)	CC 00 01 <b>21</b> 00 DD	CC 00 01 <b>21</b> 03 qb d1 d2 cs
Read Setpoint (control point)	CC 00 01 <b>70</b> 00 8E	CC 00 01 <b>70</b> 03 qb d1 d2 cs
Read Low Temperature Limit	CC 00 01 <b>40</b> 00 BE	CC 00 01 <b>40</b> 03 qb d1 d2 cs
Read High Temperature Limit	CC 00 01 <b>60</b> 00 9E	CC 00 01 <b>60</b> 03 qb d1 d2 cs
Read Cool Proportional Band (P)	CC 00 01 <b>74</b> 00 8A	CC 00 01 <b>74</b> 03 qb d1 d2 cs
Read Cool Integral (I)	CC 00 01 <b>75</b> 00 89	CC 00 01 <b>75</b> 03 qb d1 d2 cs
Read Cool Derivative (D)	CC 00 01 <b>76</b> 00 88	CC 00 01 <b>76</b> 03 qb d1 d2 cs
Read Heat Proportional Band (P)	CC 00 01 <b>71</b> 00 8D	CC 00 01 <b>71</b> 03 qb d1 d2 cs
Read Heat Integral (I)	CC 00 01 <b>72</b> 00 8C	CC 00 01 <b>72</b> 03 qb d1 d2 cs
Read Heat Derivative (D)	CC 00 01 <b>73</b> 00 8B	CC 00 01 <b>73</b> 03 qb d1 d2 cs

command bytes shown in **bold**  
 qb = qualifier byte, see Table 2  
 d1,d2 = 16 bit signed integer of the value being sent or received  
 cs = the checksum of the string (see text)  
 \* See Read Status, see Table 3

v = protocol version  
 ed = echoed command byte  
 en 01=Bad command, 03=Bad checksum  
 xx = no valid data, include in checksum

**Table 1 Commands (continued) (All bytes are in hex)**

FUNCTION	MASTER SENDS	UNIT RESPONDS
Set Setpoint (control point)*	CC 00 01 <b>F0</b> 02 d1 d2 cs	CC 00 01 <b>F0</b> 03 qb d1 d2 cs
Set Low Temp	CC 00 01 <b>C0</b> 02 d1 d2 cs	CC 00 01 <b>C0</b> 03 qb d1 d2 cs
Set High Temp	CC 00 01 <b>E0</b> 02 d1 d2 cs	CC 00 01 <b>E0</b> 03 qb d1 d2 cs
Set Cool Proportional Band (P = 0.1 - 99.9)	CC 00 01 <b>F4</b> 02 d1 d2 cs	CC 00 01 <b>F4</b> 03 qb d1 d2 cs
Set Cool Integral (I = 0 - 9.99)	CC 00 01 <b>F5</b> 02 d1 d2 cs	CC 00 01 <b>F5</b> 03 qb d1 d2 cs
Set Cool Derivative (D = 0 - 5.0)	CC 00 01 <b>F6</b> 02 d1 d2 cs	CC 00 01 <b>F6</b> 03 qb d1 d2 cs
Set Heat Proportional Band (P = 0.1 - 99.9)	CC 00 01 <b>F1</b> 02 d1 d2 cs	CC 00 01 <b>F1</b> 03 qb d1 d2 cs
Set Heat Integral (I = 0 - 9.99)	CC 00 01 <b>F2</b> 02 d1 d2 cs	CC 00 01 <b>F2</b> 03 qb d1 d2 cs
Set Heat Derivative (D = 0 - 5.0)	CC 00 01 <b>F3</b> 02 d1 d2 cs	CC 00 01 <b>F3</b> 03 qb d1 d2 cs

**Set On/Off Array**

The Set On/Off Array command is used to set and/or request unit settings. Depending on the software version, the array consists of up to 5 data bytes, d1 - d5. Each byte represents the state of one setting, see \*\* below. The master can send the array with a value of 0 (turn off/disable), 1 (turn on/enable), or 2 (do not change). 0 or 1 causes the unit to change the setting, 2 only requests the state of a setting. The slave returns values of 0 (off/disabled) or 1 (on/enabled).

(SW Version 02950.1A-E)	CC 00 01 <b>81</b> 01 d1 cs	CC 00 01 <b>81</b> 01 d1 cs
(SW Version 02950.1F-H)	CC 00 01 <b>81</b> 04 d1... d4 cs	CC 00 01 <b>81</b> 04 d1...d4 cs
(SW Version 02950.1J-K)	CC 00 01 <b>81</b> 05 d1... d5 cs	CC 00 01 <b>81</b> 05 d1...d5 cs
(SW Version 02950.1L- )	Any of the above formats	
Sample - turn unit on (02950.1J-K)	CC 00 01 <b>81</b> 05 01 02 02 02 02 6F	

\* = limited to the range of the chiller command bytes shown in **bold**  
 qb = qualifier byte, see Table 2  
 d1,d2 = 16 bit signed integer of the value being sent or received  
 cs = the checksum of the string (see text)

\*\* d1 - unit off = 0, on =1  
 d2 - external sensor disabled = 0, enabled = 1  
 d3 - fault mode disabled = 0, enabled = 1  
 d4 - tenths display disabled = 0, enabled = 1  
 d5 - Alarms on internal = 0, or external sensor =1,  
 (1 is invalid when d2=0, external sensor disabled)

**Table 2**

**QUALIFIER BYTE**

00	0 precision, no units of measure
01	0 precision, °C
10	1 precision, no units of measure
11	1 precision, °C
20	2 precision, no units of measure

Example: The temperature of 45.6°C would be represented by the qualifier 11 hex, followed by the 2 bytes 01 C8 hex (456 decimal).

**Table 3**

**READ STATUS**

BIT	d1	d2 <sup>t</sup>
b.7 = 1	External Temp Sensor Enabled	Reserved
b.6 = 1	Level 1 Warning/Auto Refill	Reserved
b.5 = 1	Low Flow Warning	Internal Temp Sensor Fault
b.4 = 1	Low Level 2 Warning	External Temp Sensor Fault
b.3 = 1	High or Low Temp Warning	High Temp/Fixed High Temp Fault
b.2 = 1	High or Low Temp Bypass	Low Temp/Fixed Low Temp Fault
b.1 = 1	Unit Faulted*	Low Flow Fault
b.0 = 1	Unit Running*	Low Level 2 Fault

t d2 is the detail of fault d1-b.1

\* mutually exclusive

## WARRANTY

Thermo Electron Corporation warrants for 24 months from date of shipment any Thermo unit according to the following terms.

Any part of the unit manufactured or supplied by Thermo and found in the reasonable judgment of Thermo to be defective in material or workmanship will be repaired at an authorized Thermo Repair Depot without charge for parts or labor. The unit, including any defective part must be returned to an authorized Thermo Repair Depot within the warranty period. The expense of returning the unit to the authorized Thermo Repair Depot for warranty service will be paid for by the buyer. Thermo's responsibility in respect to warranty claims is limited to performing the required repairs or replacements, and no claim of breach of warranty shall be cause for cancellation or rescision of the contract of sales of any unit. With respect to units that qualify for field service repairs, Thermo's responsibility is limited to the component parts necessary for the repair and the labor that is required on site to perform the repair. Any travel labor or mileage charges are the financial responsibility of the buyer.

The buyer shall be responsible for any evaluation or warranty service call (including labor charges) if no defects are found with the Thermo product.

This warranty does not cover any unit that has been subject to misuse, neglect, or accident. This warranty does not apply to any damage to the unit that is the result of improper installation or maintenance, or to any unit that has been operated or maintained in any way contrary to the operating or maintenance instructions specified in Thermo's Instruction and Operation Manual. This warranty does not cover any unit that has been altered or modified so as to change its intended use.

In addition, this warranty does not extend to repairs made by the use of parts, accessories, or fluids which are either incompatible with the unit or adversely affect its operation, performance, or durability.

Thermo reserves the right to change or improve the design of any unit without assuming any obligation to modify any unit previously manufactured.

THE FOREGOING EXPRESS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO WARRANTIES OR MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

Thermo's OBLIGATION UNDER THIS WARRANTY IS STRICTLY AND EXCLUSIVELY LIMITED TO THE REPAIR OR REPLACEMENT OF DEFECTIVE COMPONENT PARTS AND Thermo DOES NOT ASSUME OR AUTHORIZE ANYONE TO ASSUME FOR IT ANY OTHER OBLIGATION.

Thermo ASSUMES NO RESPONSIBILITY FOR INCIDENTAL, CONSEQUENTIAL, OR OTHER DAMAGES INCLUDING, BUT NOT LIMITED TO LOSS OR DAMAGE TO PROPERTY, LOSS OF PROFITS OR REVENUE, LOSS OF THE UNIT, LOSS OF TIME, OR INCONVENIENCE.

This warranty applies to units sold in the United States. Any units sold elsewhere are warranted by the affiliated marketing company of Thermo. This warranty and all matters arising pursuant to it shall be governed by the law of the State of New Hampshire, United States. All legal actions brought in relation hereto shall be filed in the appropriate state or federal courts in New Hampshire, unless waived by Thermo.